

Dashboard

Patient Search

Filter Patients

Attending Physician: All Treatment Progress: All Filter By Date: [] Filter

All Patients My Patient List

Legend

01 May 2018 - 07 May 2018

Page 1 / 2

Showing 1 - 10 of 20 10 records per page

Patients	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07	Actions
Bailey, Wayne 07 December 1965 Physician: Jean Grey	✓	!	✓	--	!!	--	✓	📷 ⋮
Diaz, Susie 07 December 1965 Physician: Charles Xavier	!	--	✓	✓	--	✓	!!	📷 ⋮
Elliott, Priscilla 07 December 1965 Physician: Jean Grey	--	✓	!	✓	!!	✓	--	📷 ⋮
Gordon, Guy 07 December 1965 Physician: Charles Xavier	✓	!	✓	--	!!	!!	!!	📷 ⋮
Mclaughlin, Chad 07 December 1965 Physician: Lemmy Kilmesiter	✓	✓	--	!!	--	✓	!	📷 ⋮
Patton, Shaun 07 December 1965 Physician: Scott Summers	✓	!	✓	--	!!	--	✓	📷 ⋮
Pope, Eddie 07 December 1965 Physician: Lemmy Kilmesiter	✓	✓	✓	✓	✓	✓	✓	📷 ⋮
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter	✓	!	✓	--	!!	--	✓	📷 ⋮
Smith, Donna 07 December 1965 Physician: Scott Summers	✓	✓	!	!!	--	✓	--	📷 ⋮
Wayne, Bruce 07 December 1965 Physician: Oswald Cobblepot	✓	!	✓	--	!!	--	✓	📷 ⋮

Page 1 / 2

Showing 1 - 10 of 20 10 records per page

Dashboard Tool Tips

Patient Search

Filter Patients

Attending Physician: All Treatment Progress: All Filter By Date: [] Filter

All Patients My Patient List

Legend

01 May 2018 - 07 May 2018

Page 1 / 2

Showing 1 - 10 of 20 10 records per page

Patients	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07	Snapshot
Bailey, Wayne 07 December 1965 Physician: Jean Grey				--		--		
Diaz, Susie 07 December 1965 Physician: Charles Xavier		--			--			
Elliott, Priscilla 07 December 1965 Physician: Jean Grey	--						--	
Gordon, Guy 07 December 1965 Physician: Charles Xavier				--				
Mclaughlin, Chad 07 December 1965 Physician: Lemmy Kilmesiter			--		--			
Patton, Shaun 07 December 1965 Physician: Scott Summers				--		--		
Pope, Eddie 07 December 1965 Physician: Lemmy Kilmesiter								
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter				--		--		
Smith, Donna 07 December 1965 Physician: Scott Summers					--		--	
Wayne, Bruce 07 December 1965 Physician: Oswald Cobblepot				--		--		

Page 1 / 2

Showing 1 - 10 of 20 10 records per page

Dashboard Filters

[Patient Search](#)

Filter Patients

Attending Physician:
 Treatment Progress:
 Filter By Date:

All P:
 Jean Grey
 Charles Xavier
 Hugo Strange
 Stephen Strange
 Bruce Banner

High Priority Flag
 Priority Flag
 Treatment Completed
 No Treatment Data
 Treatments Not Reviewed
 No Communication

May 2018

Su	Mo	Tu	We	Th	Fr	Sa
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

Page 1 / 2

10 records per page

Patients	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07	Actions
Bailey, Wayne 07 December 1965 Physician: Jean Grey	✓	! (Yellow)	✓	--	!! (Red)	--	✓	
Diaz, Susie 07 December 1965 Physician: Charles Xavier	! (Yellow)	--	✓	✓	--	✓	!! (Red)	
Elliott, Priscilla 07 December 1965 Physician: Jean Grey	--	✓	! (Yellow)	✓	!! (Red)	✓	--	
Gordon, Guy 07 December 1965 Physician: Charles Xavier	✓	! (Yellow)	✓	--	!! (Red)	!! (Red)	!! (Red)	
Mclaughlin, Chad 07 December 1965 Physician: Lemmy Kilmesiter	✓	✓	--	!! (Red)	--	✓	! (Yellow)	
Patton, Shaun 07 December 1965 Physician: Scott Summers	✓	! (Yellow)	✓	--	!! (Red)	--	✓	
Pope, Eddie 07 December 1965 Physician: Lemmy Kilmesiter	✓	✓	✓	✓	✓	✓	✓	
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter	✓	! (Yellow)	✓	--	!! (Red)	--	✓	
Smith, Donna 07 December 1965 Physician: Scott Summers	✓	✓	! (Yellow)	!! (Red)	--	✓	--	
Wayne, Bruce 07 December 1965 Physician: Oswald Cobblepot	✓	! (Yellow)	✓	--	!! (Red)	--	✓	

Page 1 / 2

Showing 1 - 10 of 20 10 records per page

Dashboard Filtered

[Patient Search](#) 🔍

Filter Patients

Attending Physician: Charles Xavier ▼
 Treatment Progress: Treatments Not Reviewed ▼
 Filter By Date: 04 May 2018
[Clear ALL Filters](#)

All Patients My Patient List

01 May 2018 - 07 May 2018

[Legend](#)

Page 1 / 1

Showing 1 - 5 of 5 10 records per page

Patients	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07	Actions
Diaz, Susie 07 December 1965 Physician: Charles Xavier	!	--	✓	✓	--	✓	!!	
Gordon, Guy 07 December 1965 Physician: Charles Xavier	✓	!	✓	--	!!	!!	!!	

Page 1 / 1

Showing 1 - 5 of 5 10 records per page

Dashboard My Patients Tab

[Patient Search](#) 🔍

Filter Patients

Attending Physician: ▼
 Treatment Progress: ▼
 Filter By Date:

◀ 01 May 2018 - 07 May 2018 ▶

[Legend](#)

◀ ◁ Page 1 / 1 ▷ ▶

Showing 1 - 3 of 3 records per page

Patients	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07	Actions
Bailey, Wayne 07 December 1965 Physician: Jean Grey	✓	! (Yellow)	✓	--	!! (Red)	--	✓	
Mclaughlin, Chad 07 December 1965 Physician: Lemmy Kilmesiter	✓	✓	--	!! (Red)	--	✓	! (Yellow)	
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter	✓	! (Yellow)	✓	--	!! (Red)	--	✓	

◀ ◁ Page 1 / 1 ▷ ▶

Showing 1 - 3 of 3 records per page

Dashboard Pagination

Patient Search

Filter Patients

Attending Physician: All
 Treatment Progress: All
 Filter By Date:

All Patients My Patient List

Legend

01 May 2018 - 07 May 2018

Page 1 / 2 Showing 1 - 10 of 20 records per page

Patients		Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07		
Bailey, Wayne 07 December 1965 Physician: Jean Grey	-	✓	!	✓	--	!!	--	✓	📷	☰
Diaz, Susie 07 December 1965 Physician: Charles Xavier	+	!	--	✓	✓	--	✓	!!	📷	☰
Elliott, Priscilla 07 December 1965 Physician: Jean Grey	+	--	✓	!	✓	!!	✓	--	📷	☰
Gordon, Guy 07 December 1965 Physician: Charles Xavier	+	✓	!	✓	--	!!	!!	!!	📷	☰
Mclaughlin, Chad 07 December 1965 Physician: Lemmy Kilmesiter	-	✓	✓	--	!!	--	✓	!	📷	☰
Patton, Shaun 07 December 1965 Physician: Scott Summers	+	✓	!	✓	--	!!	--	✓	📷	☰
Pope, Eddie 07 December 1965 Physician: Lemmy Kilmesiter	+	✓	✓	✓	✓	✓	✓	✓	📷	☰
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter	-	✓	!	✓	--	!!	--	✓	📷	☰
Smith, Donna 07 December 1965 Physician: Scott Summers	+	✓	✓	!	!!	--	✓	--	📷	☰
Wayne, Bruce 07 December 1965 Physician: Oswald Cobblepot	+	✓	!	✓	--	!!	--	✓	📷	☰

Page 1 / 2 Showing 1 - 10 of 20 records per page

Dashboard Legend

[Patient Search](#)

Filter Patients

Attending Physician:
 Treatment Progress:
 Filter By Date:

All Patients My Patient List

[Legend](#)

01 May 2018 - 07 May 2018

Legend

Dashboard Notifications

- Treatment Completed
- High Priority Flag
- No Treatment Data
- Multiple (3+) Treatments
- Treatment Reviewed
- Priority Flag
- No Communication

Dashboard Actions


- Patient Snapshot
- Add to My Patients List
- Device Settings
- Remove From My Patients List

Showing 1 - 10 of 20 records per page

	Wed 04	Thurs 05	Fri 06	Sat 07	Actions	
07 December 1965 Physician: Lemmy Kilmesiter						
Reed, Adrian 07 December 1965 Physician: Lemmy Kilmesiter						
Smith, Donna 07 December 1965 Physician: Scott Summers						
Wayne, Bruce 07 December 1965 Physician: Oswald Cobblepot						

Patient Search

Patient Actions

 [Return to Dashboard](#)














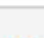






Patient Search

Patient Name

Patient Status All ▾ Search

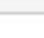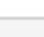

 Page 1 ▾ / 2  
Showing 1 - 10 of 20 10 ▾ records per page

Patient Name		Birthdate	Baxter ID	Physician	Actions
Bailey, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Benjamin, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	 
Bradley, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bruce, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	 
Bailey, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Benjamin, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	 
Bradley, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bruce, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	 
Bradley, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bruce, Wayne		07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	 



 Page 1 ▾ / 2  
Showing 1 - 10 of 20 10 ▾ records per page

Patient Search Program New Device

Patient Actions

Return to Dashboard

Patient Search

Patient Name

Patient Status All

Page 1 / 2 10 records per page

Patient Name	Birthdate	Device ID	Dr. Name	Actions
Bailey, Wayne	07 Decem		Dr. Barron, Isabell	Program New Device
Benjamin, Wayne	07 Decem		Dr. Barron, Isabell	Program New Device
Bradley, Wayne	07 Decem		Dr. Barron, Isabell	Program New Device
Bruce, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bailey, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Benjamin, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bradley, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bruce, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bradley, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device
Bruce, Wayne	07 December 1965	O-VAL01-1000079	Dr. Barron, Isabell	Program New Device

Program New Device

Select Device to program

Amia

Home Choice Claria

Page 1 / 2 Showing 1 - 10 of 20 10 records per page

Snapshot 7 Days

- Patient Actions**
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

[<< Previous Patient](#)
[Next Patient >>](#)
Smith, John

Attending Physician: Dr. Barron, Isabell

Date of Birth: 22 November 1977

Clinic Patient ID:

Baxter Patient ID: O-TRAIN-1003042





Homechoice Claria

 Select Timeframe 7 Days
01 May 2018 - 07 May 2018
[Legend](#)

	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07
				--		--	

Treatment Data	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07
Program Name	ROUTINE	GREEN/RED	ROUTINE	ROUTINE	ROUTINE	--	ROUTINE
Total UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Manual UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Day UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Night UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Initial Drain Volume (mL)	33	13	3	76	52	--	51
Pre Weight (kg)	130.7	129.8	130.7	129.8	130.7	--	129.8
Post Weight (kg)	132.7	129.5	125.6	127.1	130.0	--	128.1
Pre Blood Pressure (mmHg)	136/70	126/56	119/61	119/60	132/71	--	131/66
Pre Blood Glucose	226 mg/dl	236 mg/dl	225 mg/dl	219 mg/dl	221 mg/dl	--	226 mg/dl
Post Blood Glucose	237 mg/dl	238 mg/dl	233 mg/dl	226 mg/dl	240 mg/dl	--	250 mg/dl

Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.

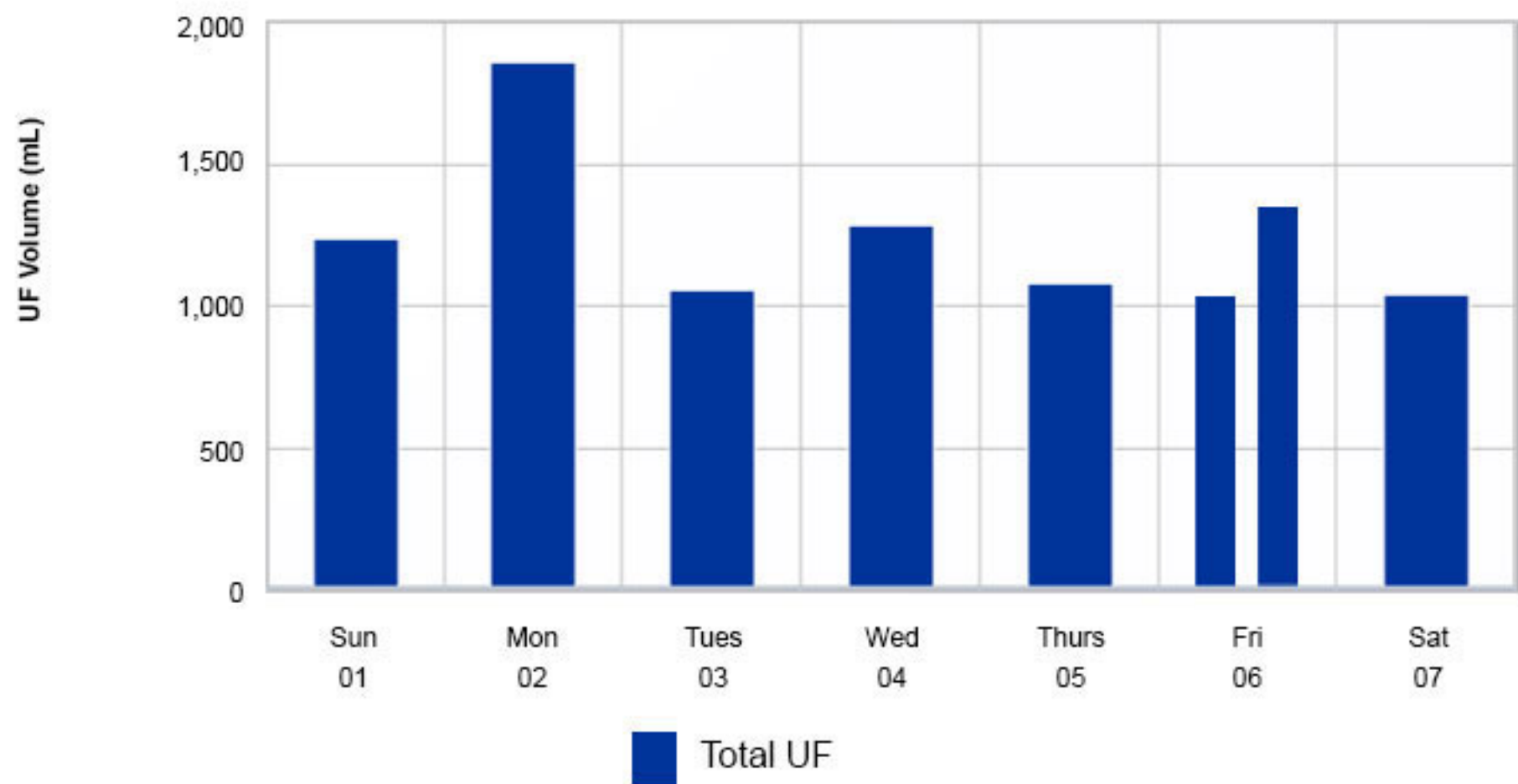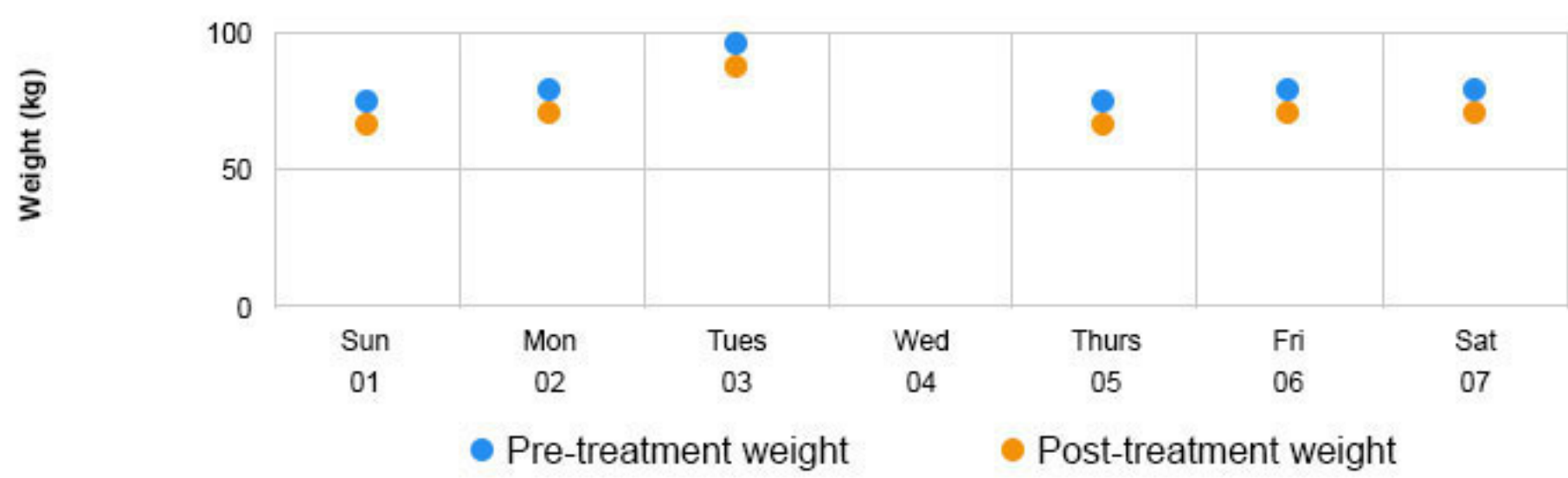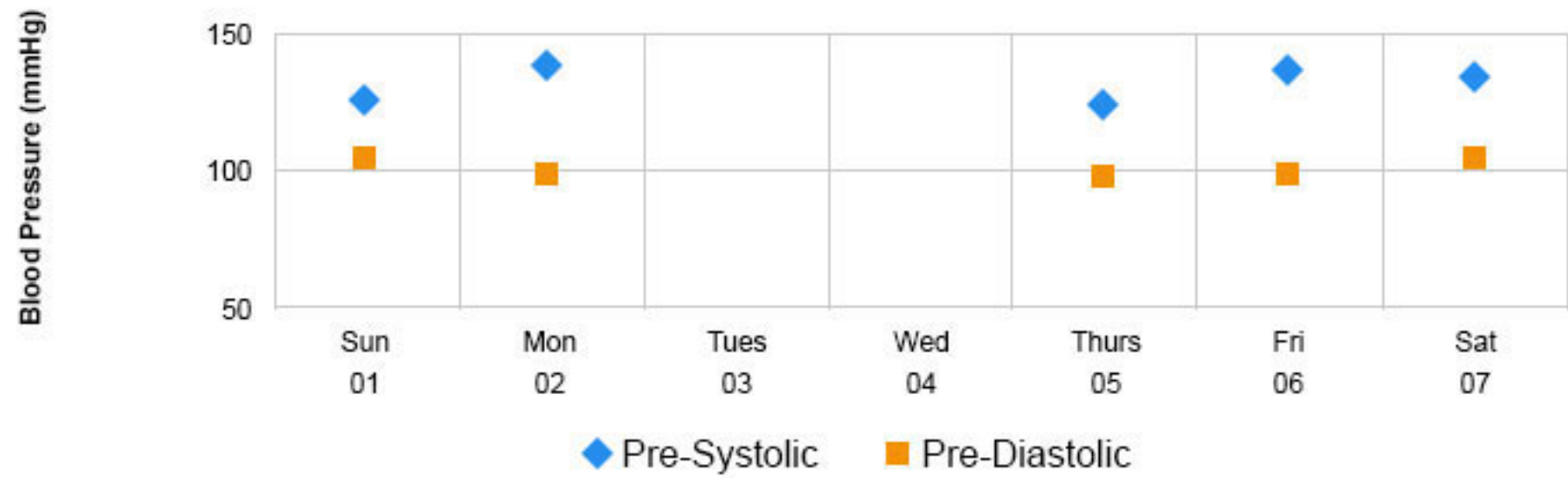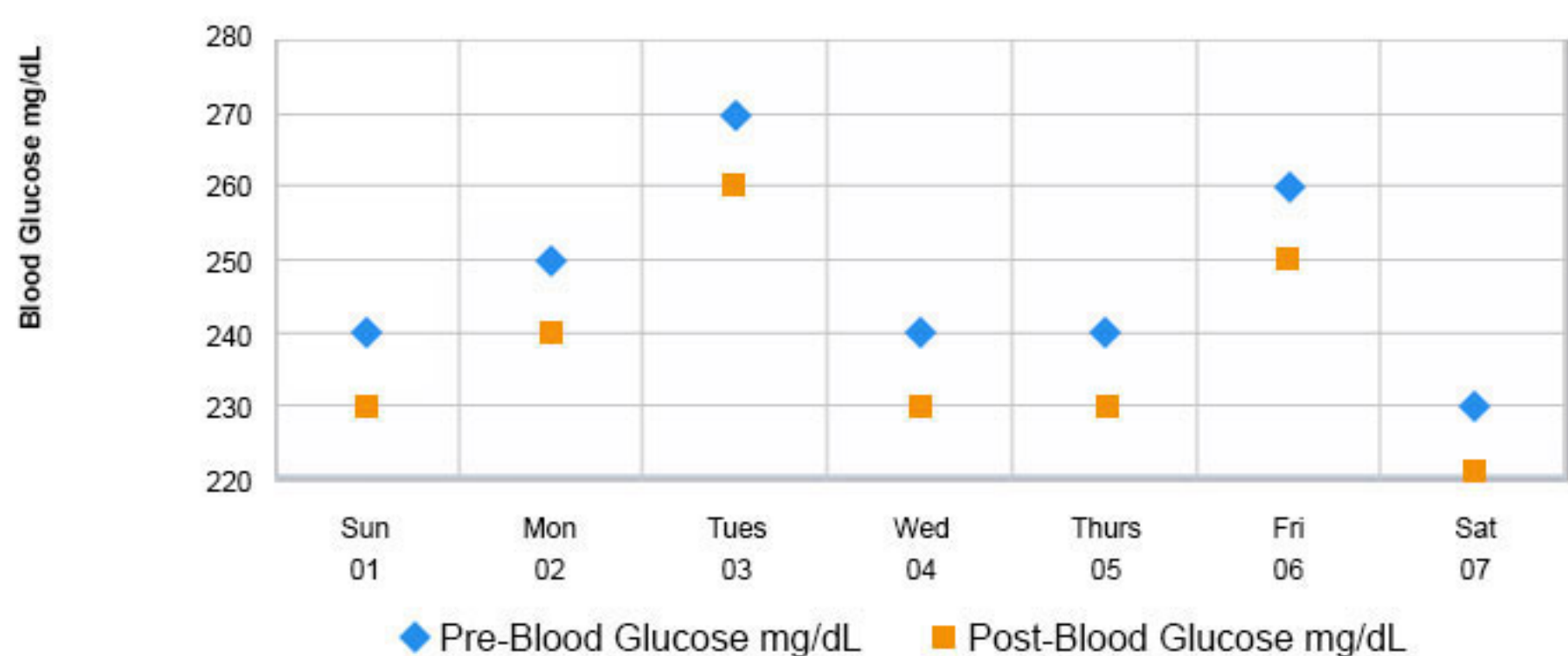Night Cycle UF

Weight

Pre Blood Pressure

Blood Glucose


Snapshot 7 Days Rollover & Flyout

- Patient Actions**
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria

Select Timeframe 7 Days





01 May 2018 - 07 May 2018

[Legend](#)

	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07
	✔	!	✔	--	!!	--	✔

Treatment Data	Sun 01	Mon 02	Tues 03	Wed 04	Thurs 05	Fri 06	Sat 07
Program Name	ROUTINE	GREEN/RED	ROUTINE	ROUTINE	ROUTINE	--	ROUTINE
Total UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Manual UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Day UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Night UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Initial Drain Volume (mL)	33	13	3	76	52	--	51
Pre Weight (kg)	130.7	129.8	130.7	129.8	130.7	--	129.8
Post Weight (kg)	132.7	129.5	125.6	127.1	130.0	--	128.1
Pre Blood Pressure (mmHg)	136/70	126/56	119/61	119/60	132/71	--	131/66
Pre Blood Glucose	226 mg/dl	236 mg/dl	225 mg/dl	219 mg/dl	221 mg/dl	--	226 mg/dl
Post Blood Glucose	237 mg/dl	238 mg/dl	233 mg/dl	226 mg/dl	240 mg/dl	--	250 mg/dl

! Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.

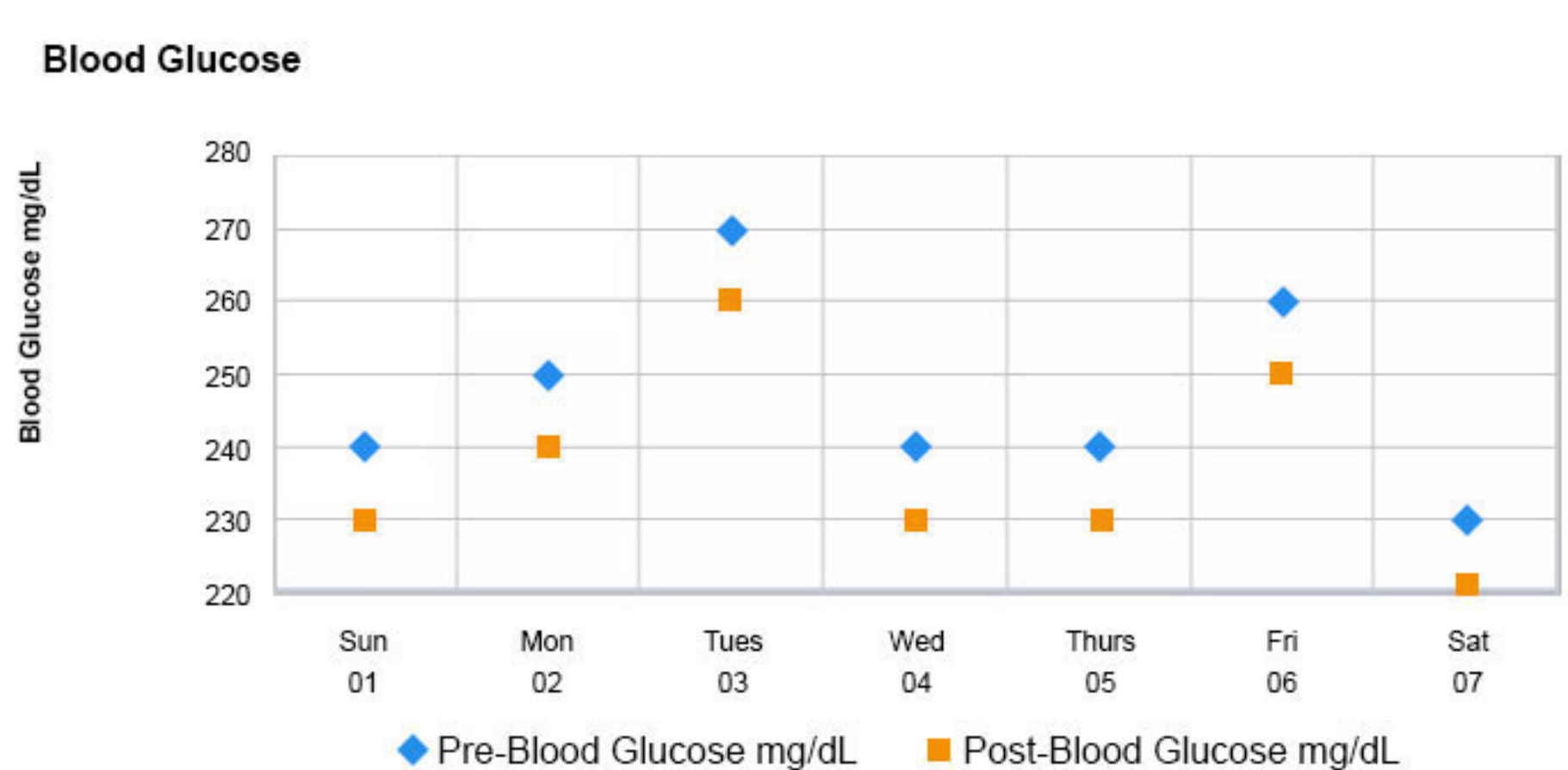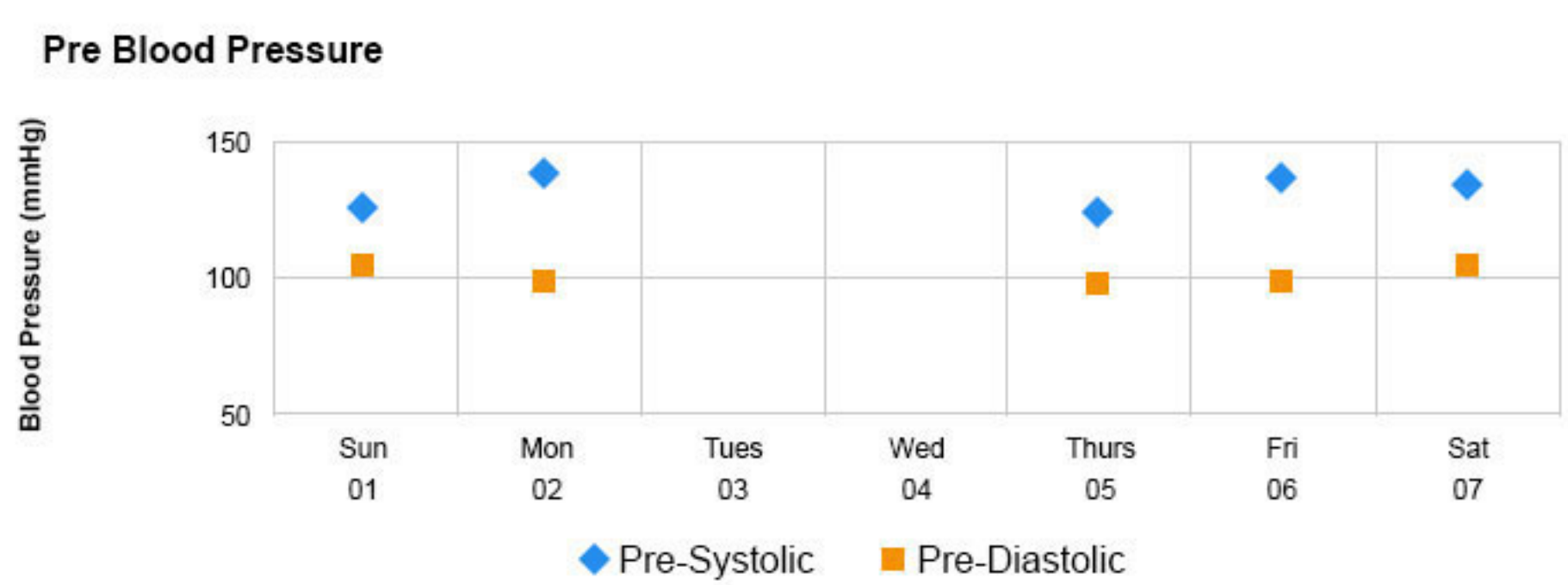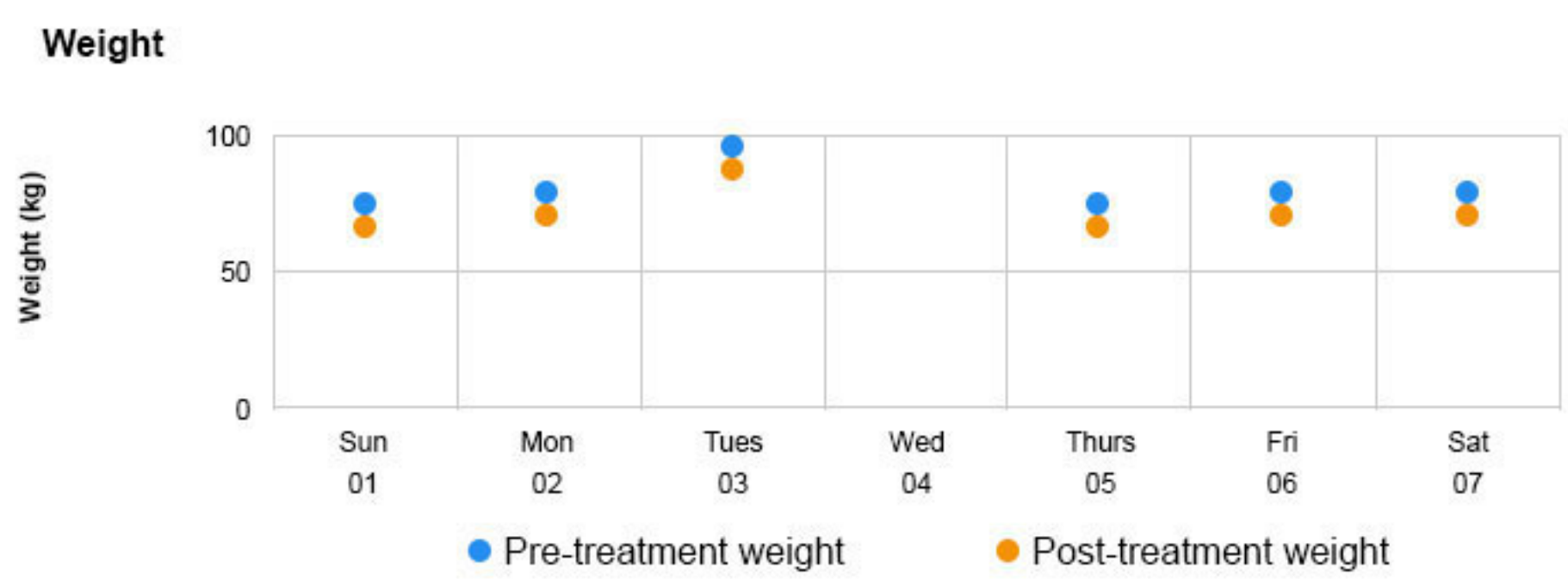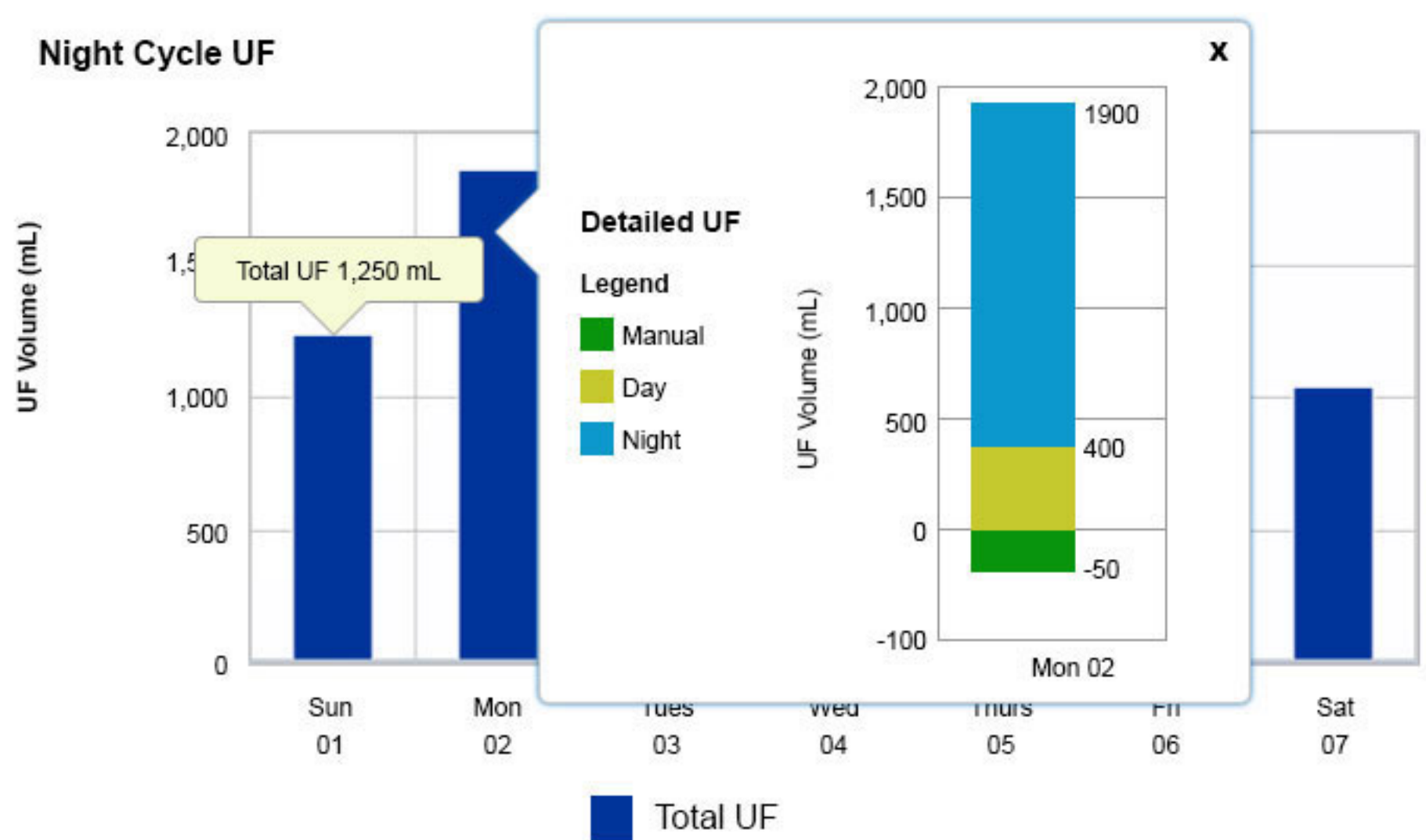

Snapshot 7 Days Legend

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient

Next Patient >>

Smith, John

Attending Physician: Dr. Barron, Isabell

Date of Birth: 22 November 1977

Clinic Patient ID:

Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria

Select Timeframe 7 Days


01 May 2018 - 07 May 2018

Legend


Treatment Data							Sat 07
Program Name							ROUTINE
Total UF (mL)							1,0333
Manual UF (mL)							1,0333
Day UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Night UF (mL)	1,231	1,857	1,054	1,279	1,073	--	1,0333
Initial Drain Volume (mL)	33	13	3	76	52	--	51
Pre Weight (kg)	130.7	129.8	130.7	129.8	130.7	--	129.8
Post Weight (kg)	132.7	129.5	125.6	127.1	130.0	--	128.1
Pre Blood Pressure (mmHg)	136/70	126/56	119/61	119/60	132/71	--	131/66
Pre Blood Glucose	226 mg/dl	236 mg/dl	225 mg/dl	219 mg/dl	221 mg/dl	--	226 mg/dl
Post Blood Glucose	237 mg/dl	238 mg/dl	233 mg/dl	226 mg/dl	240 mg/dl	--	250 mg/dl

Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.


Night Cycle UF

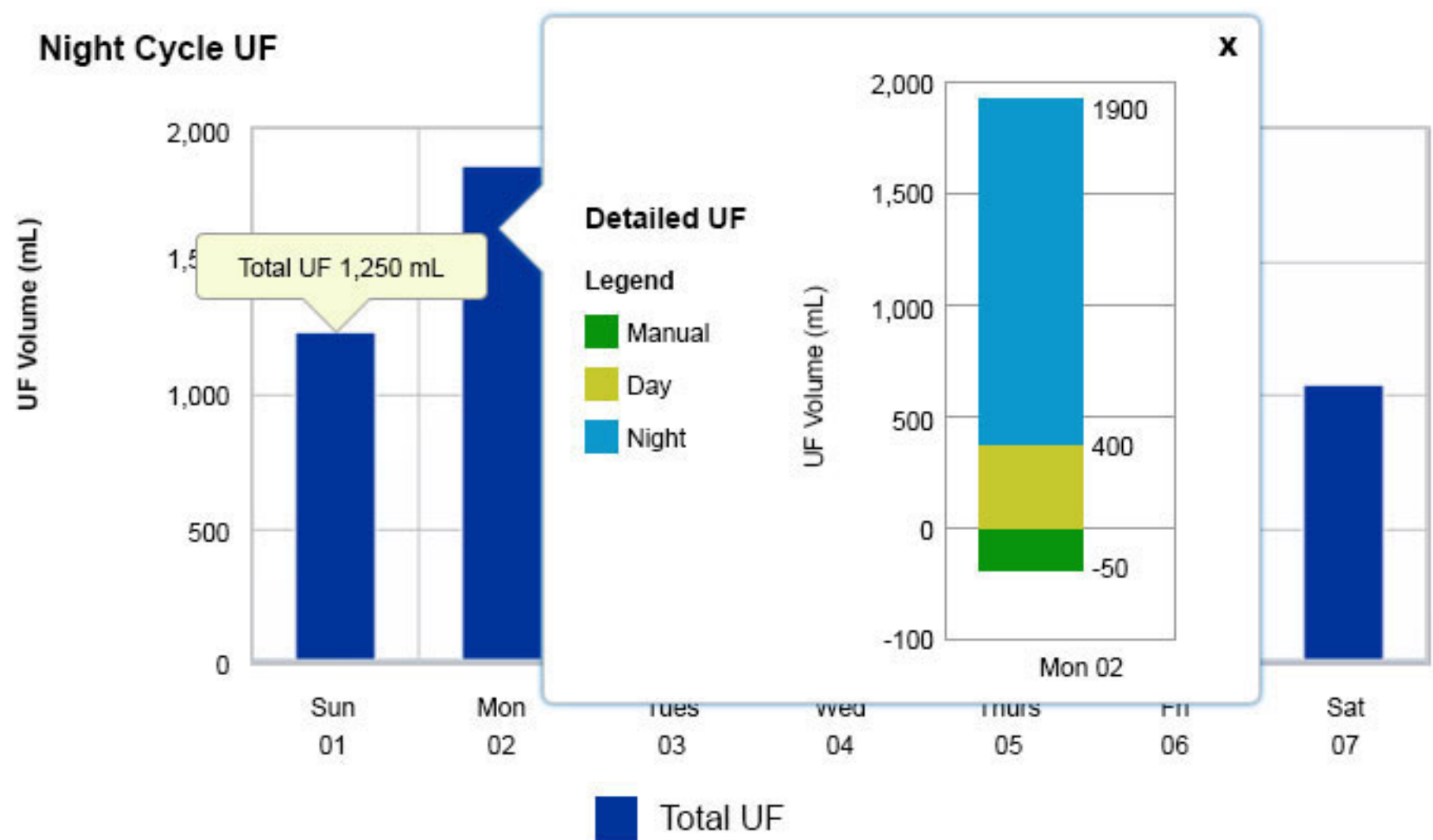


Weight

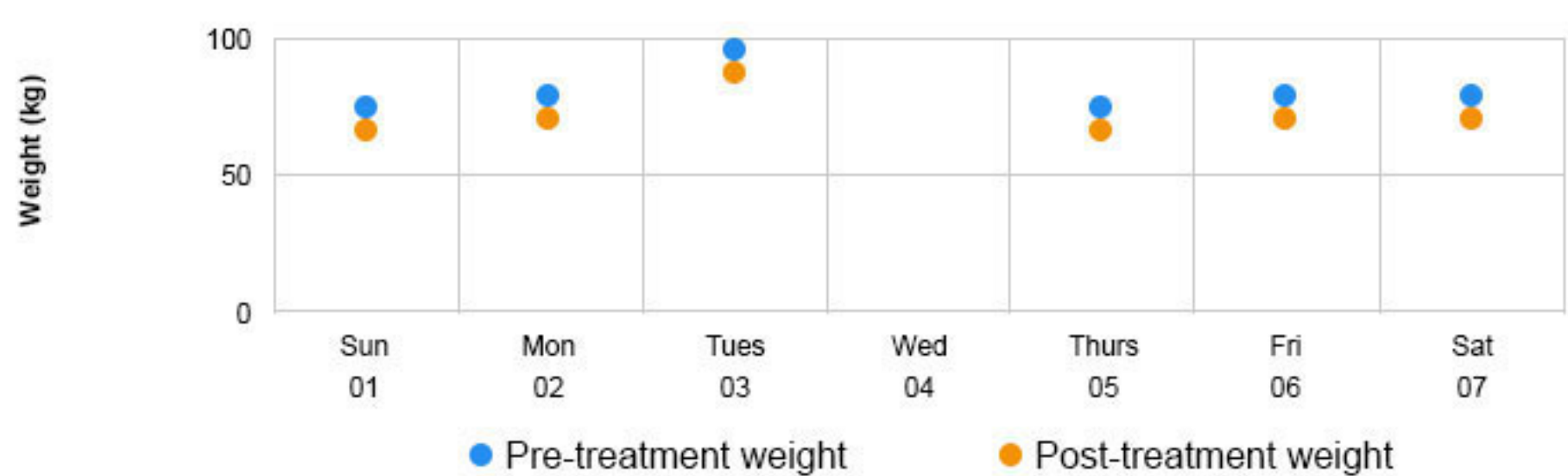

Pre Blood Pressure

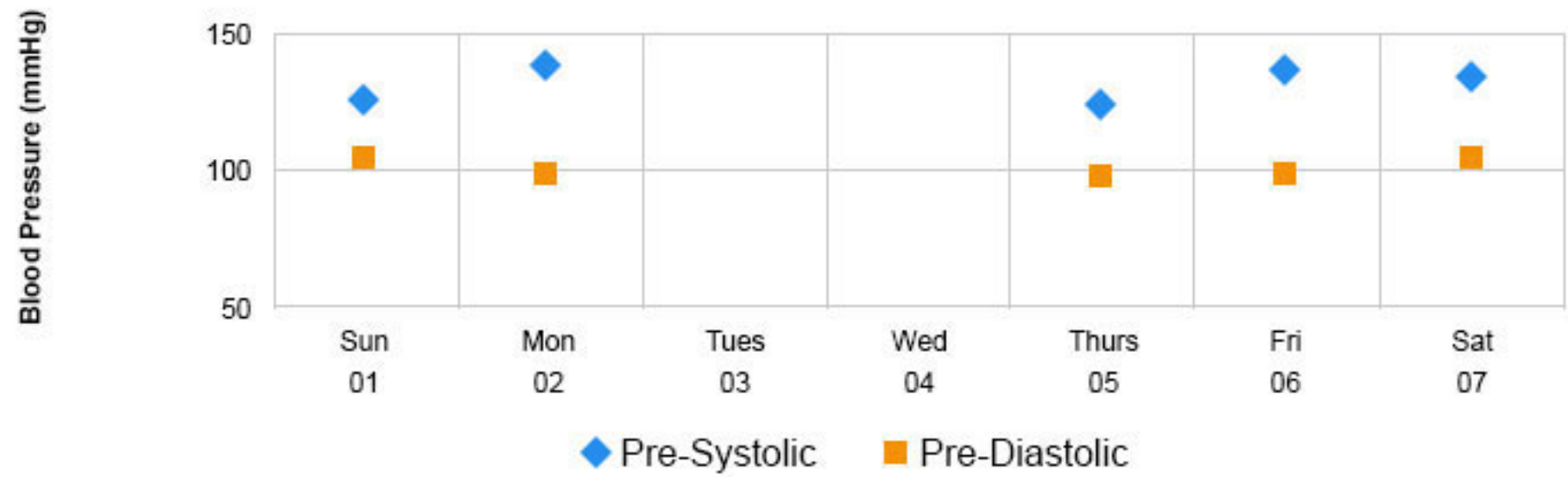

Blood Glucose

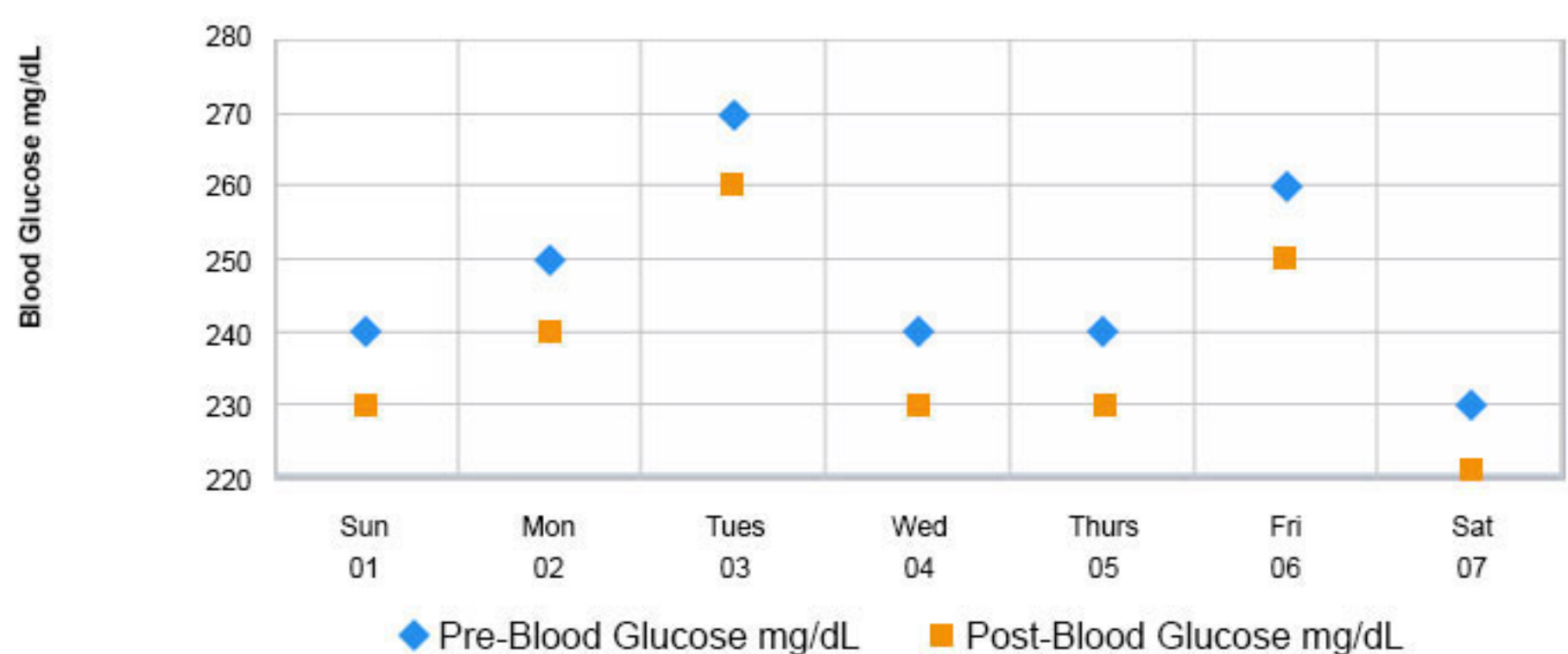

Snapshot 30 Days

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient

Next Patient >>

Smith, John

Attending Physician: Dr. Barron, Isabell

Date of Birth: 22 November 1977

Clinic Patient ID:

Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria


Select Timeframe 30 Days

01 May 2018 - 07 May 2018




! Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.


Night Cycle UF

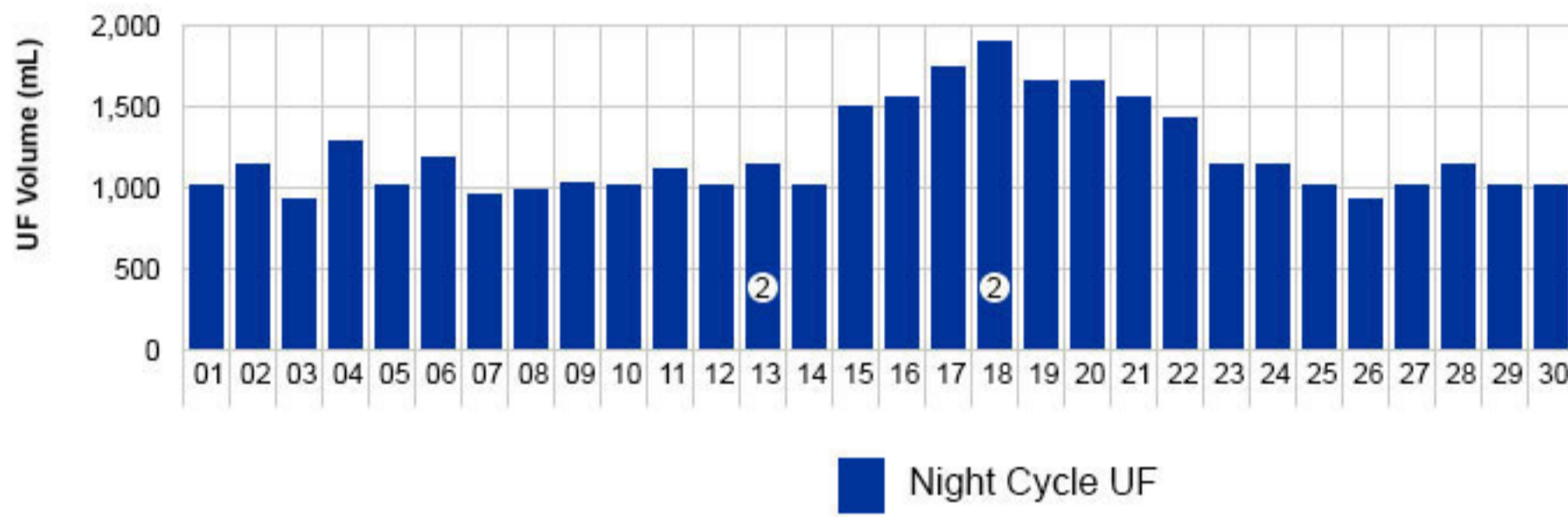


Weight

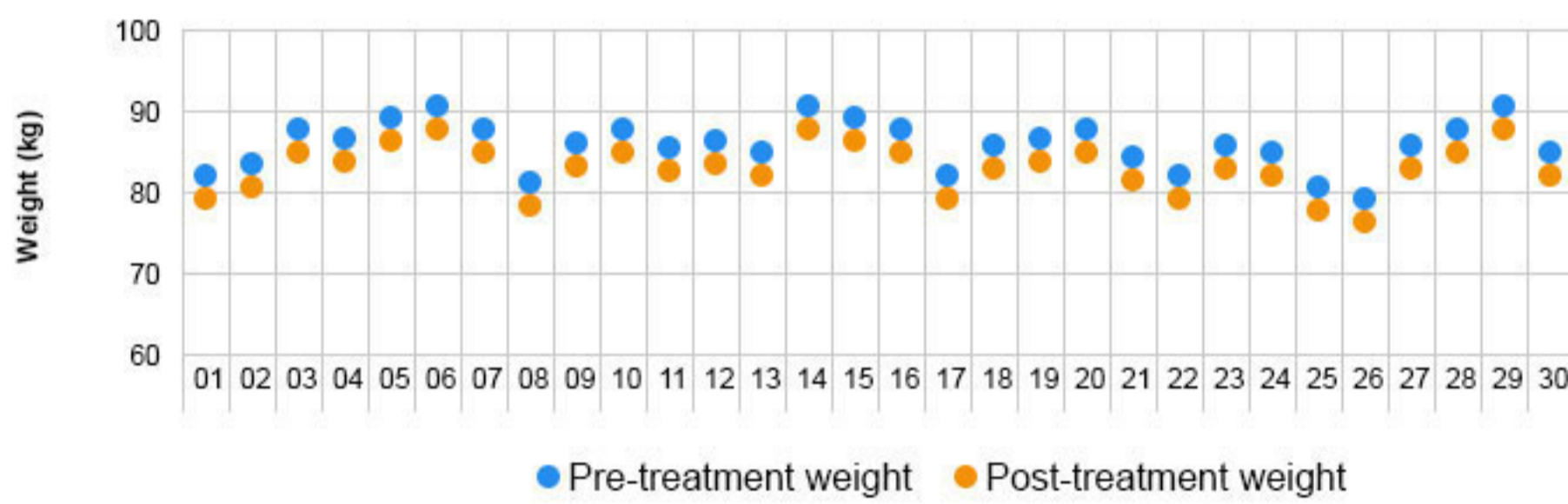

Pre-Blood Pressure (mmHg)

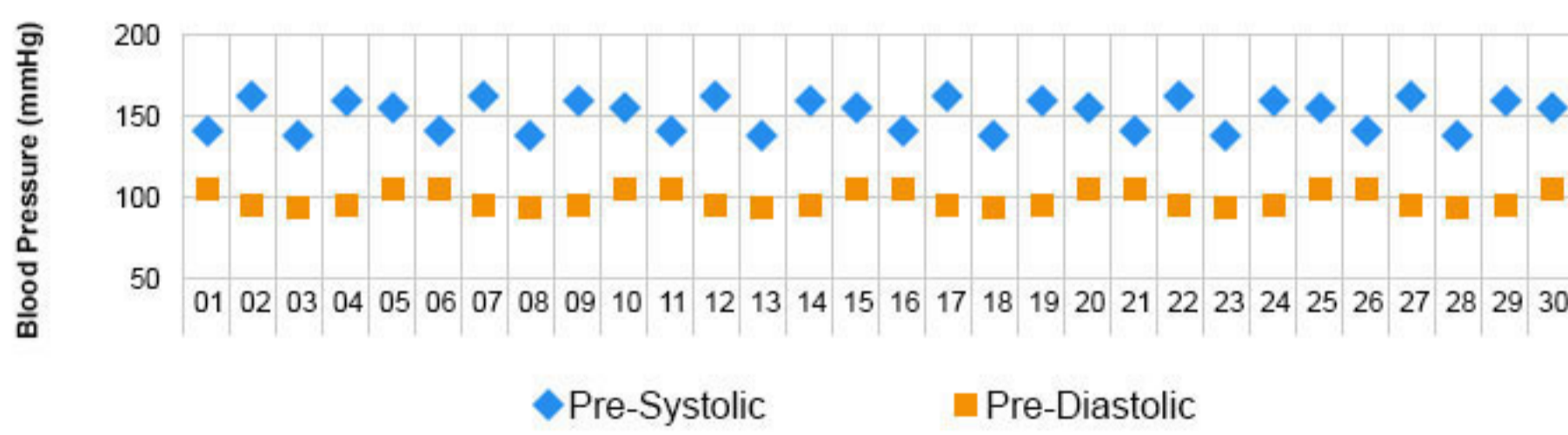

Blood Glucose

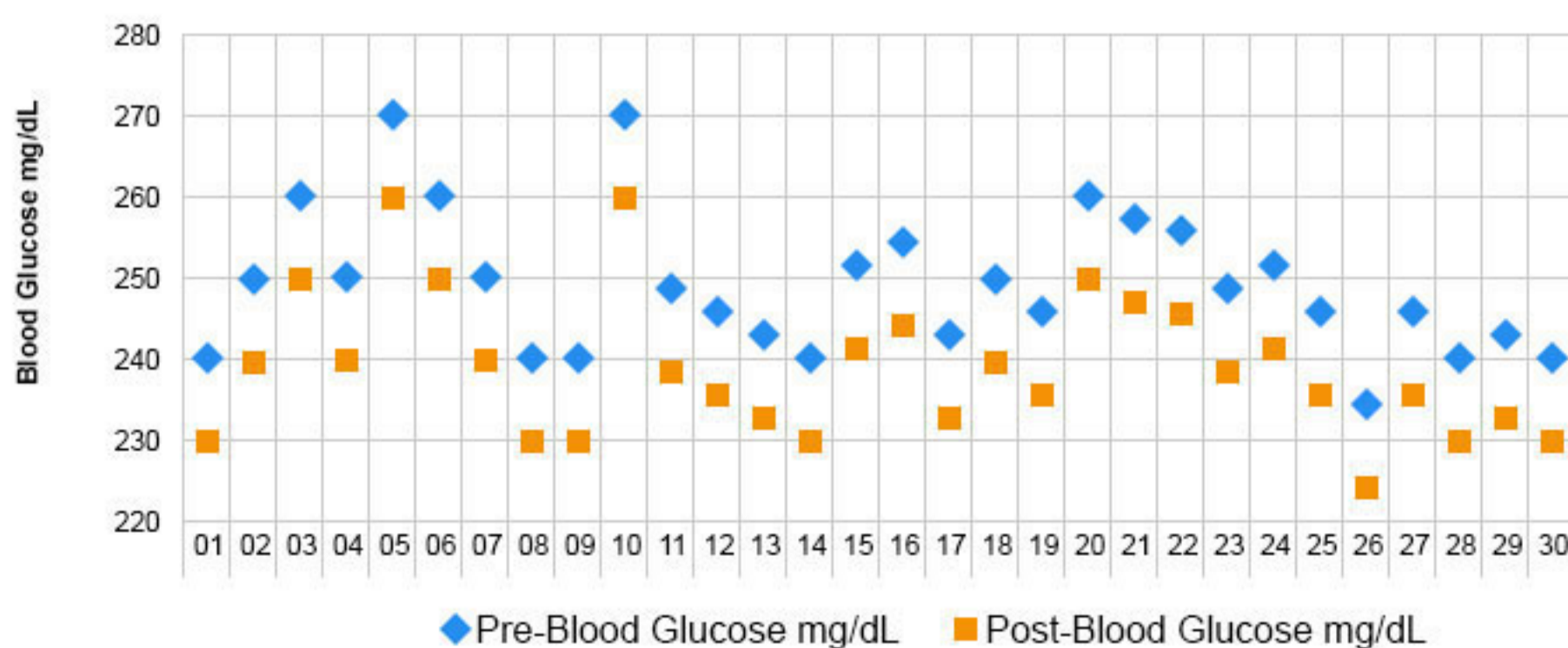

**Snapshot 30 Days
Single & Double Rollover**

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient


Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
 Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042





Homechoice Claria

Select Timeframe 30 Days

01 May 2018 - 07 May 2018

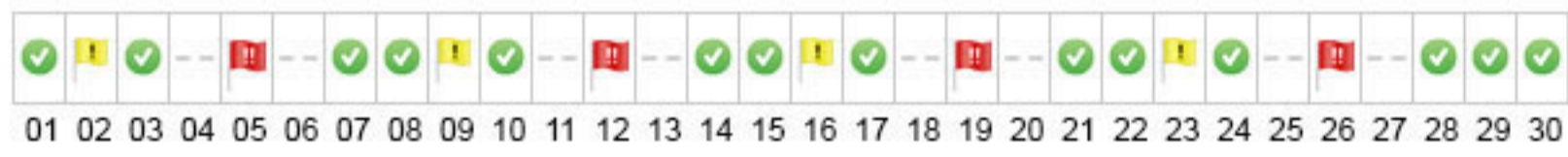

! Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.

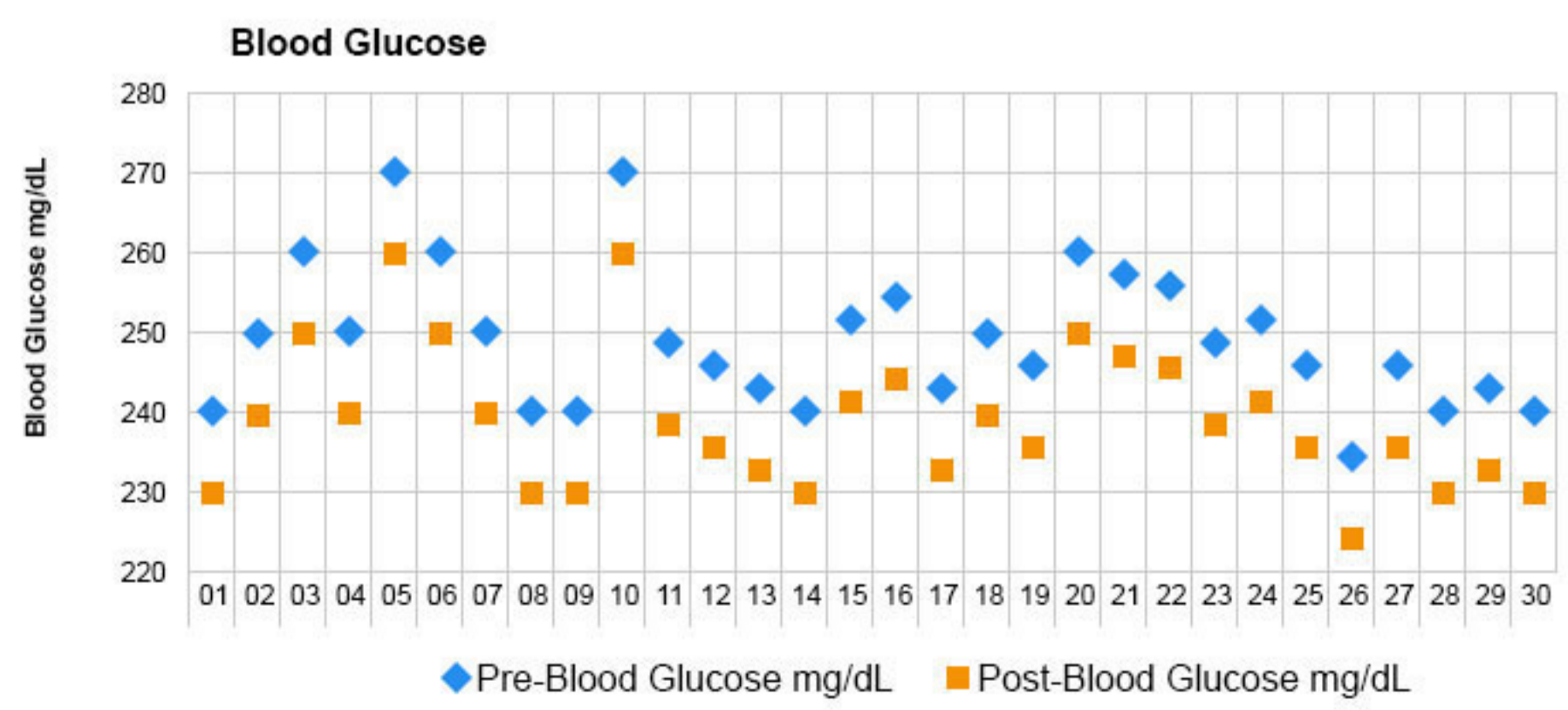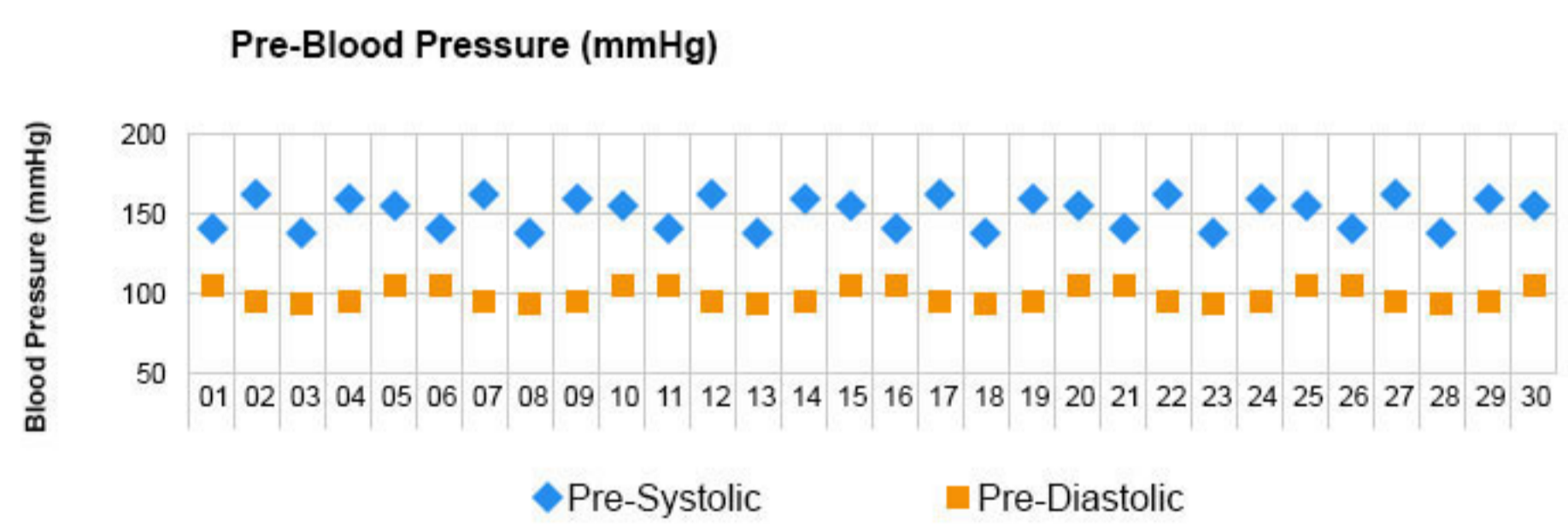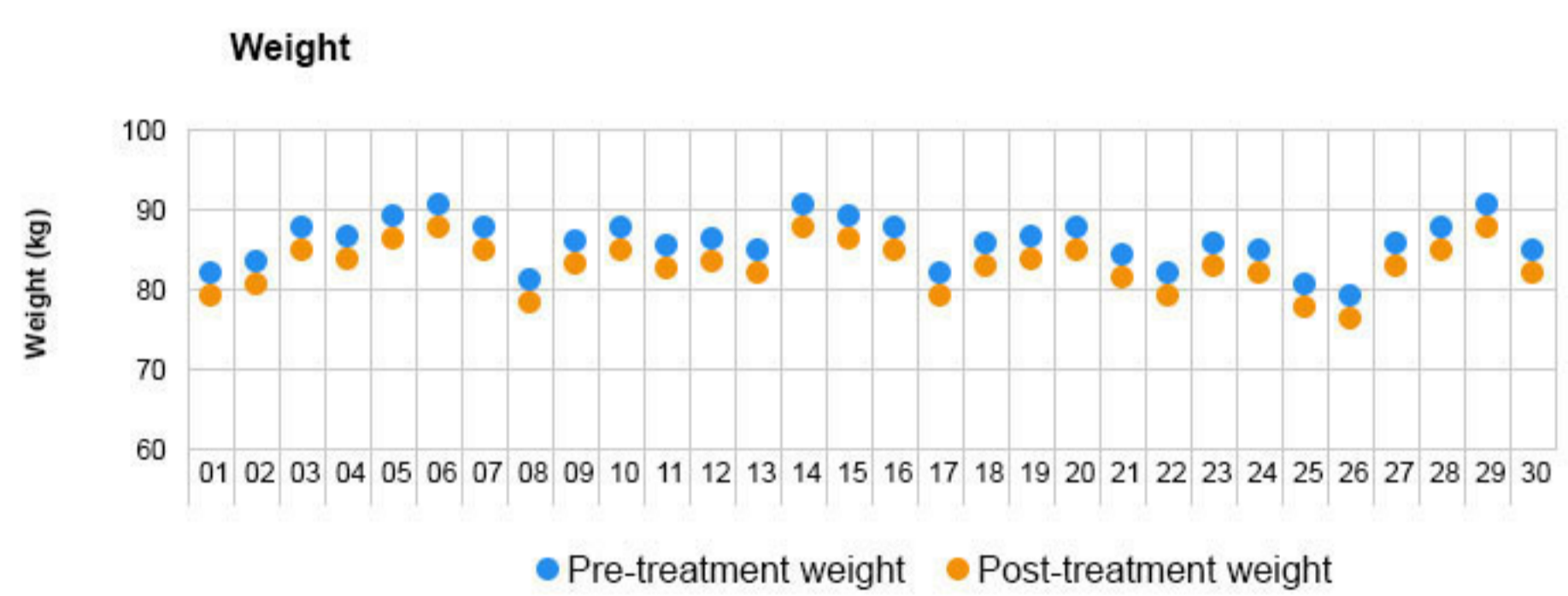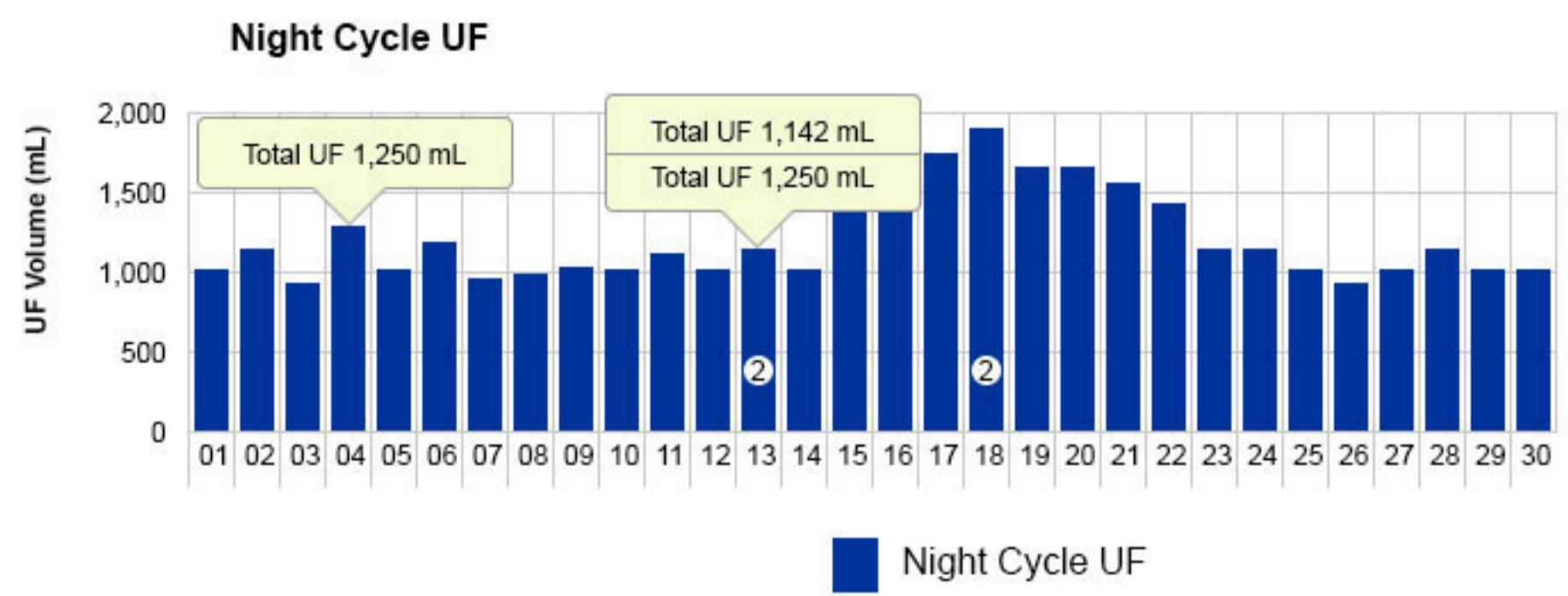

Snapshot 30 Days Single Flyout

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient


Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
 Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042




Homechoice Claria

Select Timeframe 30 Days

01 May 2018 - 07 May 2018

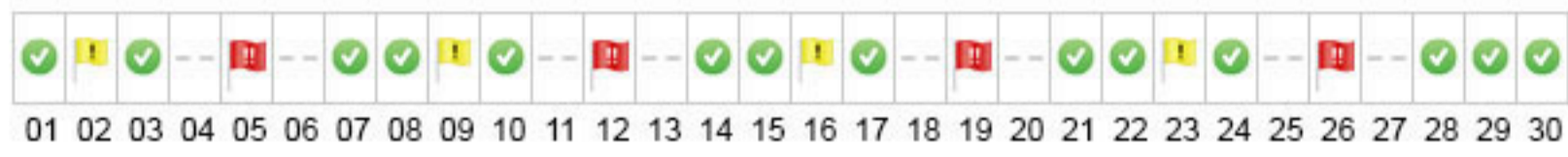

! Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.

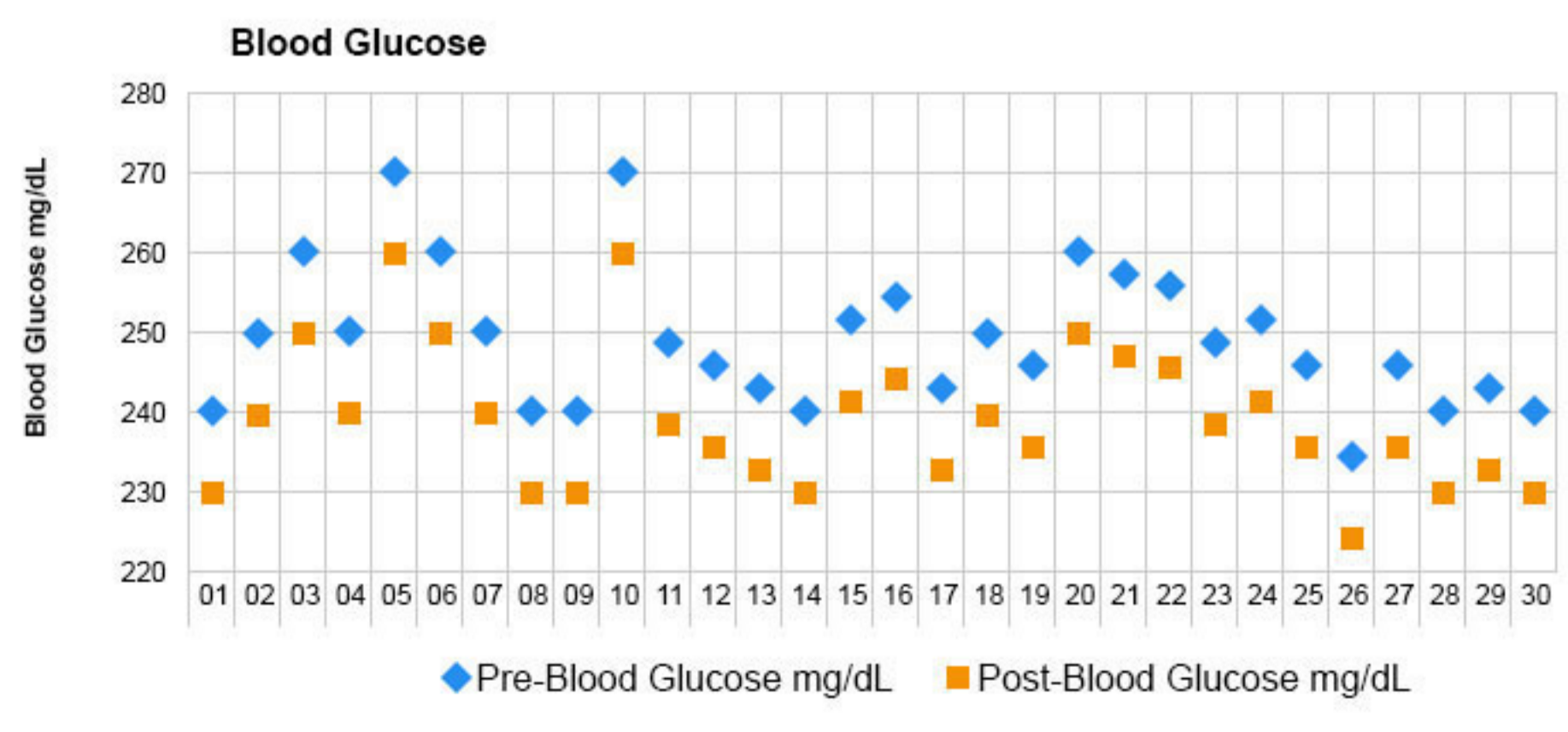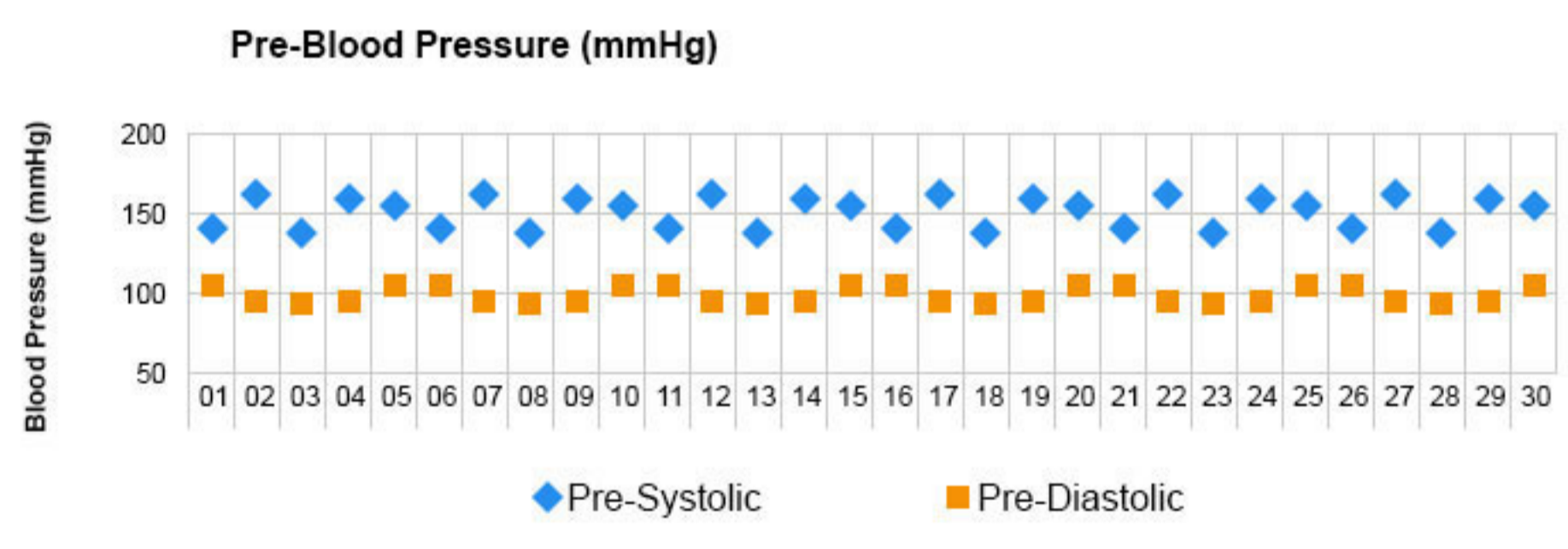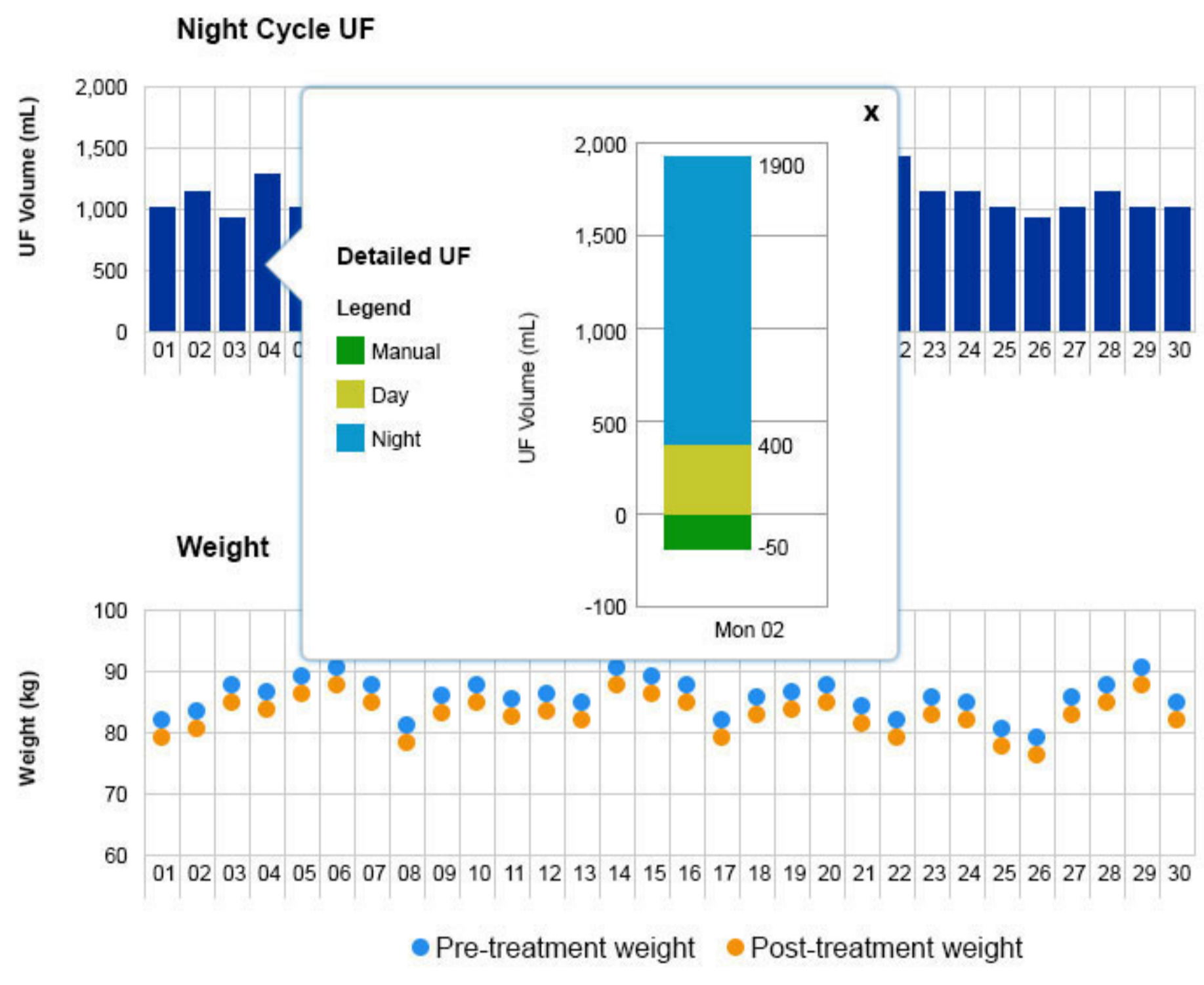

**Snapshot 30 Days
Double Flyout**

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient


Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria



Select Timeframe 30 Days

01 May 2018 - 07 May 2018

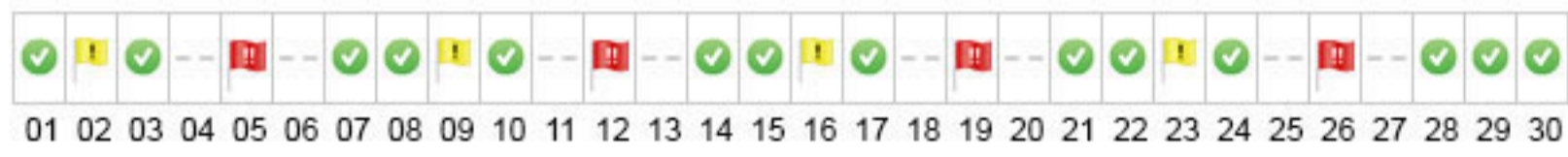


! Rollover the graphs for additional details. Click on the graph in Total Filtration for a detailed breakdown of the data.


Night Cycle UF

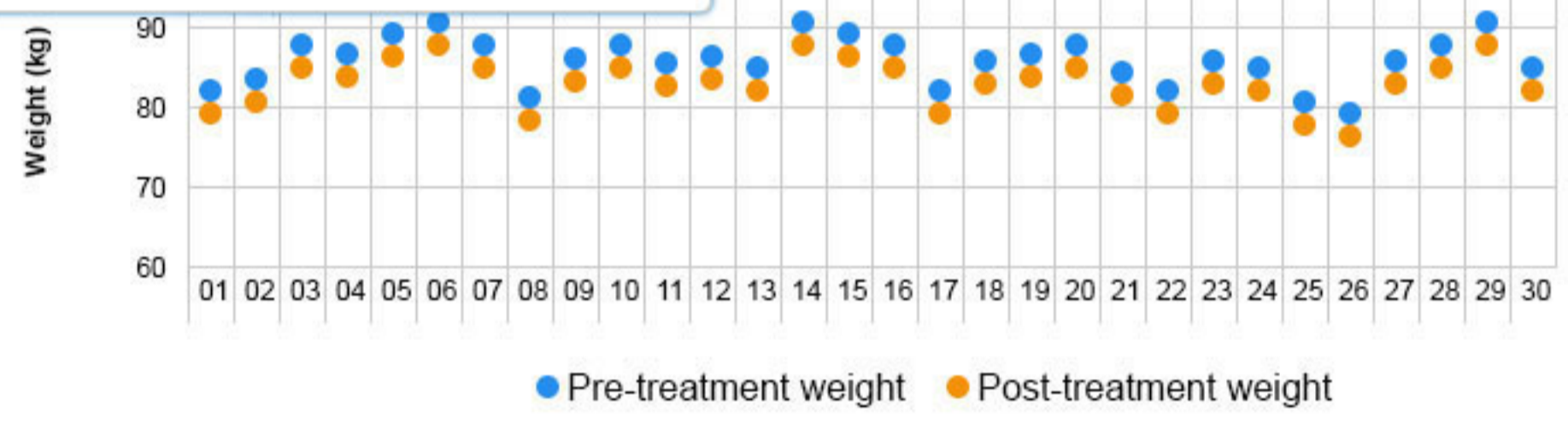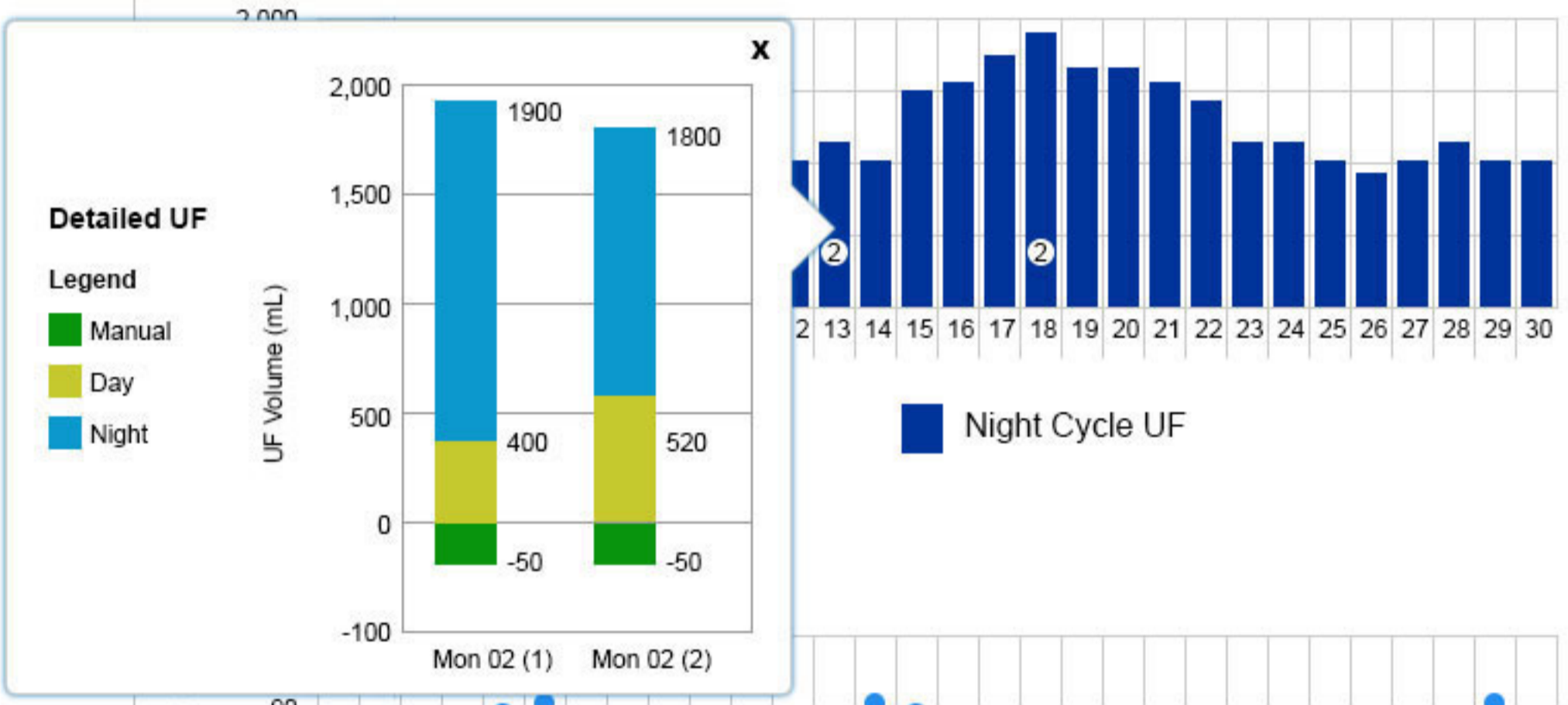

Pre-Blood Pressure (mmHg)

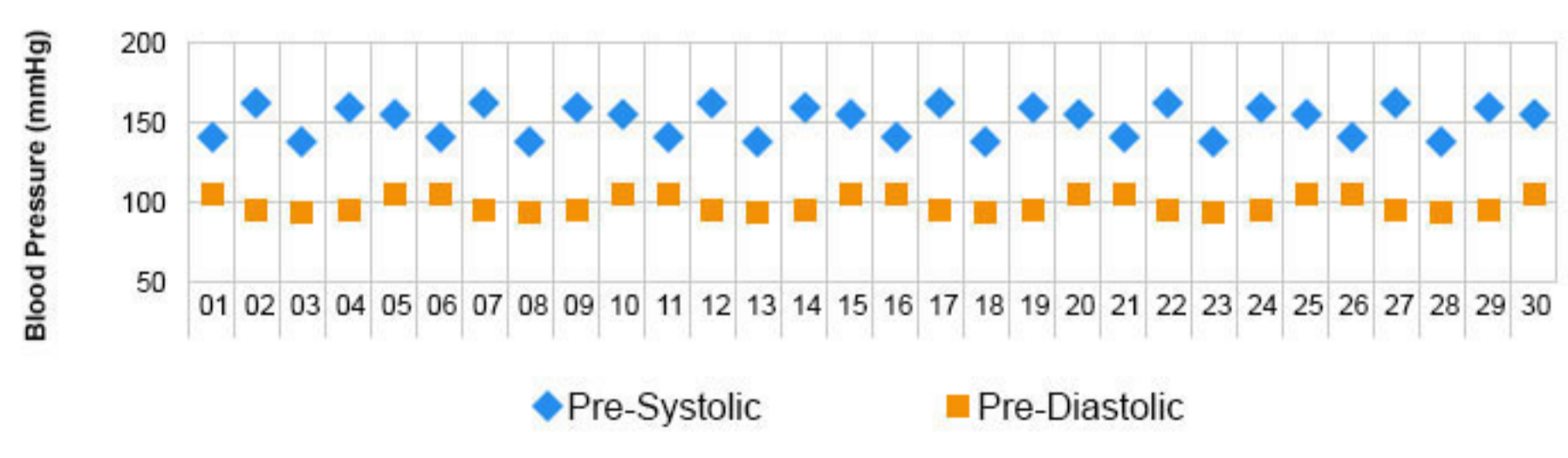

Blood Glucose

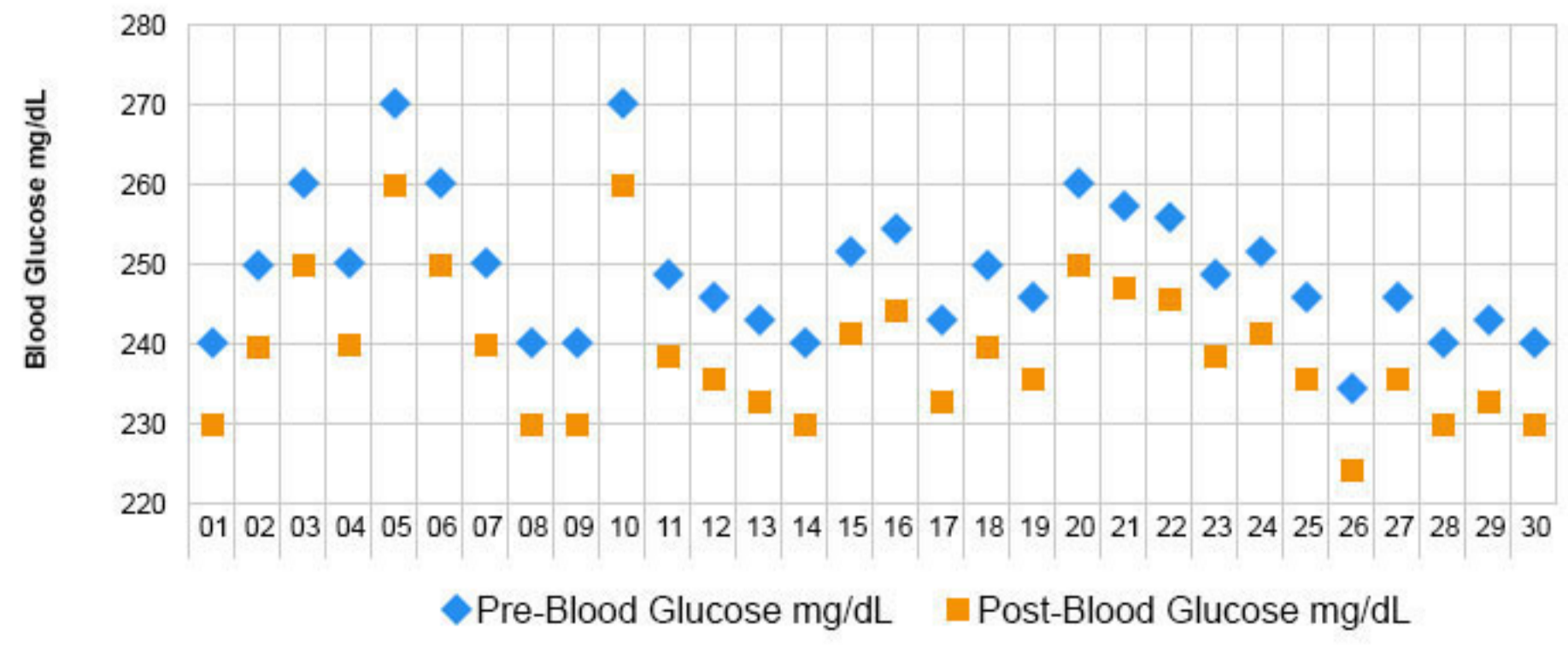

Treatment Summary

- Patient Actions**
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

[<< Previous Patient](#)
[Next Patient >>](#)

Smith, John Attending Physician: Dr. Barron, Isabell
 Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria [Generate PDF Report](#)

[<< Previous Treatment](#)
02 Apr 2018, 10:45 am
[Next Treatment >>](#)

- Fill/Dwell Bypass Count: 3
- Events During Treatment: 7


Alarms (4) and Deviations (1)

Vitals, Solutions, Total Therapy UF

Prescribed Device Program Settings – Device Program Name: HiDose							
Type/Therapy Time (HH:MM)	Number of Cycles	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
Night Time-Based / 09:00	4	2,400	--	--	--	--	--

Actual Therapy Details							
Time Stamp (HH:MM:SS)	Cycle	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
23:01:20	Initial Drain	--	--	--	0:05	52	-248
23:06:50	Night Cycle 1	2,399	0:12	1:36	0:32	2,279	-120
01:27:02	Night Cycle 2	2,400	0:13	1:19	0:36	3,044	644
03:34:34	Night Cycle 3	2,400	0:12	1:26	0:31	2,470	70
05:43:24	Night Cycle 4	2,400	0:13	1:27	0:24	2,879	479
07:46:22	Last Fill	300	0:03	--	--	--	--

Cycler Total							
08:41	--	9,899	0:53	5:48	2:08	10,724	w

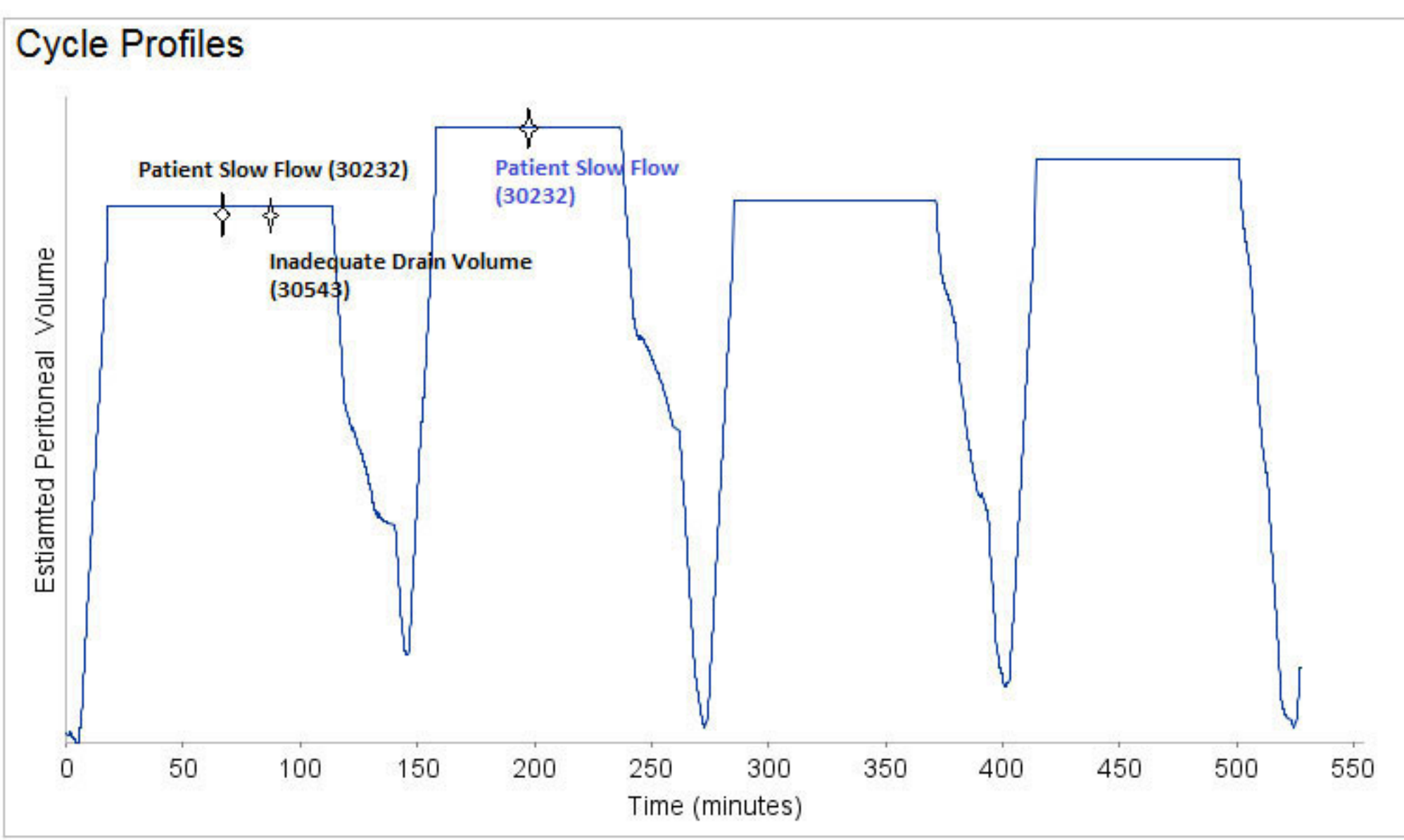

Manual Exchange	Solution Type	Fill Volume	Drain Volume	UF
1	Dianeal	ml	ml	ml
2	Dianeal	ml	ml	ml

Device Program Details

Device Data, Drug Survey

To mark this treatment as reviewed click the Treatment Reviewed button on the right. [Treatment Reviewed](#)

Treatment Summary Everything Expanded

- Patient Actions**
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient | Next Patient >>

Smith, John | Date of Birth: 22 November 1977 | Clinic Patient ID: | Attending Physician: Dr. Barron, Isabell
 Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria | Generate PDF Report

<< Previous Treatment | 02 Apr 2018, 10:45 am | Next Treatment >>

Fill/Dwell Bypass Count: 3 | Events During Treatment: 7

Alarms (4) and Deviations (1)

Time Stamp (HH:MM:SS)	Event Description	Deviations
01:40:22	Patient Line is Missing (30643)	Lost Dwell
01:40:22	Patient Slow Flow (30232)	
01:40:22	Inadequate Drain Volume (30543)	
01:40:22	Patient Slow Flow (30232)	

Vitals, Solutions, Total Therapy UF

Actual PD Solutions			Day Therapy UF	
Red Line	--	--	Manual Exchange UF	--
White 1 Line	Dianeal Low Calcium, 2.50%	6,000 mL	Initial Drain UF	-248 mL
White 2 Line	Dianeal Low Calcium, 2.50%	6,000 mL	Recovered Initial Drain UF	--
Blue Line	--	--	Day Cycle UF	--
Solution Total	--	12,000 mL	Total Day UF	-248 mL
Vitals			Night Therapy UF	
	Pre-Treatment	Post-Treatment	Night Cycle UF	1,073 mL
Weight	130.7 kg	Not Surveyed	Recovered Last Drain UF	--
Blood Pressure	132/71 mmHg	Not Surveyed	Total Night UF	1,073 mL
Pulse	88 beats/min	Not Surveyed	Total Therapy UF	
Blood Glucose	226 mg/dL	237 mg/dL	Total UF	825 mL
Temperature	Not Surveyed	Not Surveyed		

Prescribed Device Program Settings - Device Program Name: HiDose


Type/Therapy Time (HH:MM)	Number of Cycles	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
Night Time-Based / 09:00	4	2,400	--	--	--	--	--

Actual Therapy Details

Time Stamp (HH:MM:SS)	Cycle	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
23:01:20	Initial Drain	--	--	--	0:05	52	-248
23:06:50	Night Cycle 1	2,399	0:12	1:36	0:32	2,279	-120
01:27:02	Night Cycle 2	2,400	0:13	1:19	0:36	3,044	644
03:34:34	Night Cycle 3	2,400	0:12	1:26	0:31	2,470	70
05:43:24	Night Cycle 4	2,400	0:13	1:27	0:24	2,879	479
07:46:22	Last Fill	300	0:03	--	--	--	--

Cycler Total							
08:41	--	9,899	0:53	5:48	2:08	10,724	w

Cycle Profiles

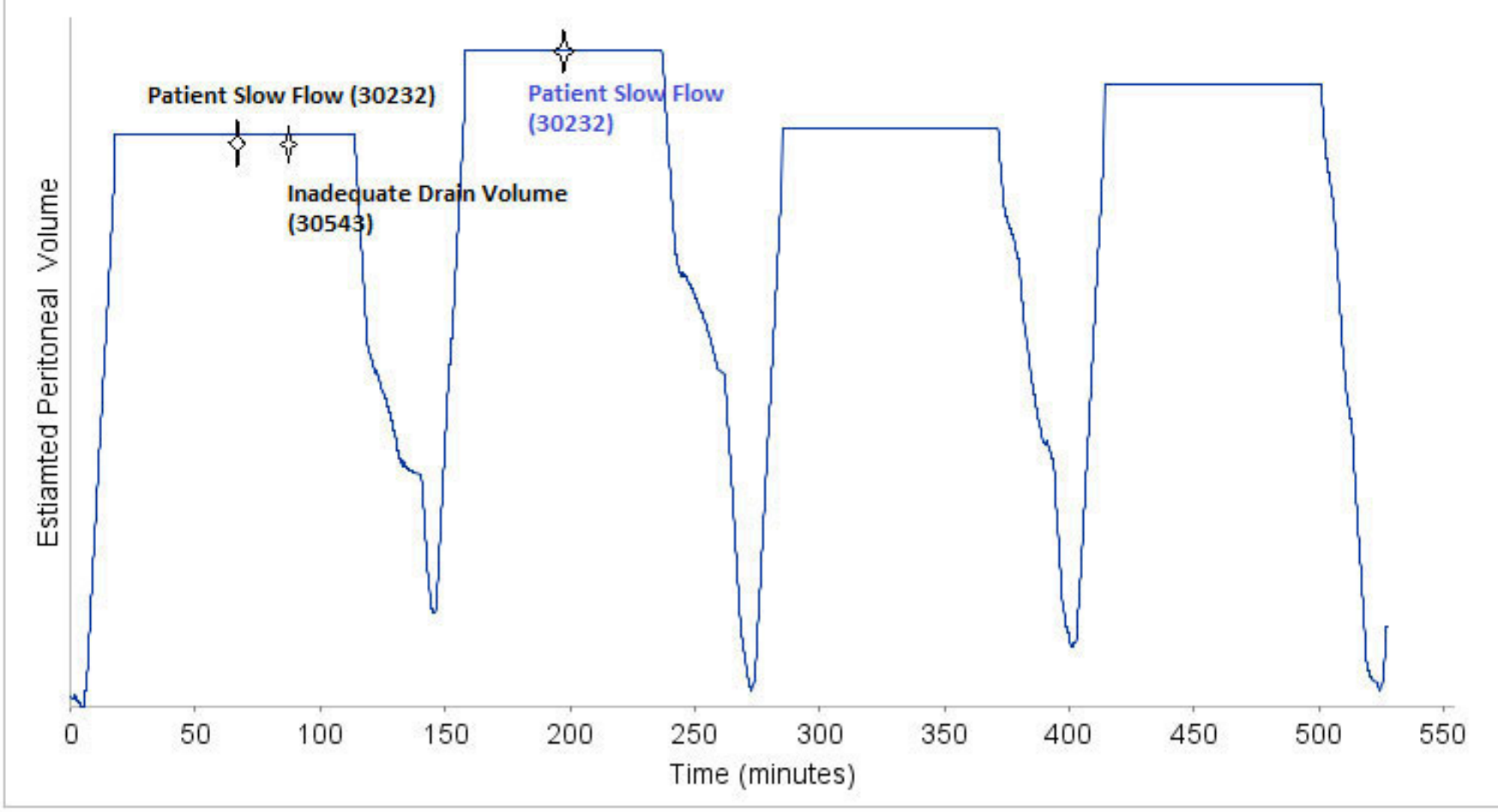

Manual Exchange	Solution Type	Fill Volume	Drain Volume	UF
1	Dianeal	ml	ml	ml
2	Dianeal	ml	ml	ml

Device Program Details

Total Therapy Volume	9,900 mL	Tidal Reserve Volume	--	Extra Last Drain Mode	Always
Last Fill Volume	300 mL	Max Peritoneal Volume	3,600 mL	Extra Drain UF Limit	--
Last Solution Type	Same	Estimated Night UF	1,100 mL	Extra Last Drain Alert	Off
Full Drain Frequency	--	Min Initial Drain Volume	0 mL	Min Initial Drain Time	5 minutes
Tidal Percent	--	Min Night Drain Time	20 minutes	Min Day Drain Volume Percent	--
Tidal Fill Volume	--	Effluent Sample Reminder	Disabled	Min Day Drain Time	--

Device Data, Drug Survey

Device Data		Drugs	Applied
Device ID	27060118	ESA	Not Surveyed
Device Software Version	2.2.0-5	Antibiotics	Not Surveyed
		Heparin	Not Surveyed
		Vitamin D	Not Surveyed

To mark this treatment as reviewed click the Treatment Reviewed button on the right. | Treatment Reviewed

Treatment Summary Calendar Flyout

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Treatment Summary
- Claria Device Settings

<< Previous Patient

Next Patient >>



Attending Physician: Dr. Barron, Isabell
Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

02 Apr 2018, 10:45 am

Generate PDF Report
Next Treatment >>

ts During Treatment: 7

Time Stamp (HH:MM:SS)	Event Description
01:40:22	Patient Line is Missing (30643)
01:40:22	Patient Slow Flow (30232)
01:40:22	Inadequate Drain Volume (30543)
01:40:22	Patient Slow Flow (30232)

Deviations
Lost Dwell

Vitals, Solutions, Total Therapy UF

Actual PD Solutions		
Red Line	--	--
White 1 Line	Dianeal Low Calcium, 2.50%	6,000 mL
White 2 Line	Dianeal Low Calcium, 2.50%	6,000 mL
Blue Line	--	--
Solution Total	--	12,000 mL

Day Therapy UF	
Manual Exchange UF	--
Initial Drain UF	-248 mL
Recovered Initial Drain UF	--
Day Cycle UF	--
Total Day UF	-248 mL

Vitals	Pre-Treatment	Post-Treatment
Weight	130.7 kg	Not Surveyed
Blood Pressure	132/71 mmHg	Not Surveyed
Pulse	88 beats/min	Not Surveyed
Blood Glucose	226 mg/dL	237 mg/dL
Temperature	Not Surveyed	Not Surveyed

Night Therapy UF	
Night Cycle UF	1,073 mL
Recovered Last Drain UF	--
Total Night UF	1,073 mL

Total Therapy UF	
Total UF	825 mL

Prescribed Device Program Settings - Device Program Name: HiDose


Type/Therapy Time (HH:MM)	Number of Cycles	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
Night Time-Based / 09:00	4	2,400	--	--	--	--	--

Actual Therapy Details

Time Stamp (HH:MM:SS)	Cycle	Fill Volume (mL)	Fill Time (h:mm)	Dwell Time (h:mm)	Drain Time (h:mm)	Drain Volume (mL)	UF/Cycle (mL)
23:01:20	Initial Drain	--	--	--	0:05	52	-248
23:06:50	Night Cycle 1	2,399	0:12	1:36	0:32	2,279	-120
01:27:02	Night Cycle 2	2,400	0:13	1:19	0:36	3,044	644
03:34:34	Night Cycle 3	2,400	0:12	1:26	0:31	2,470	70
05:43:24	Night Cycle 4	2,400	0:13	1:27	0:24	2,879	479
07:46:22	Last Fill	300	0:03	--	--	--	--

Cycler Total							
08:41	--	9,899	0:53	5:48	2:08	10,724	w

Cycle Profiles

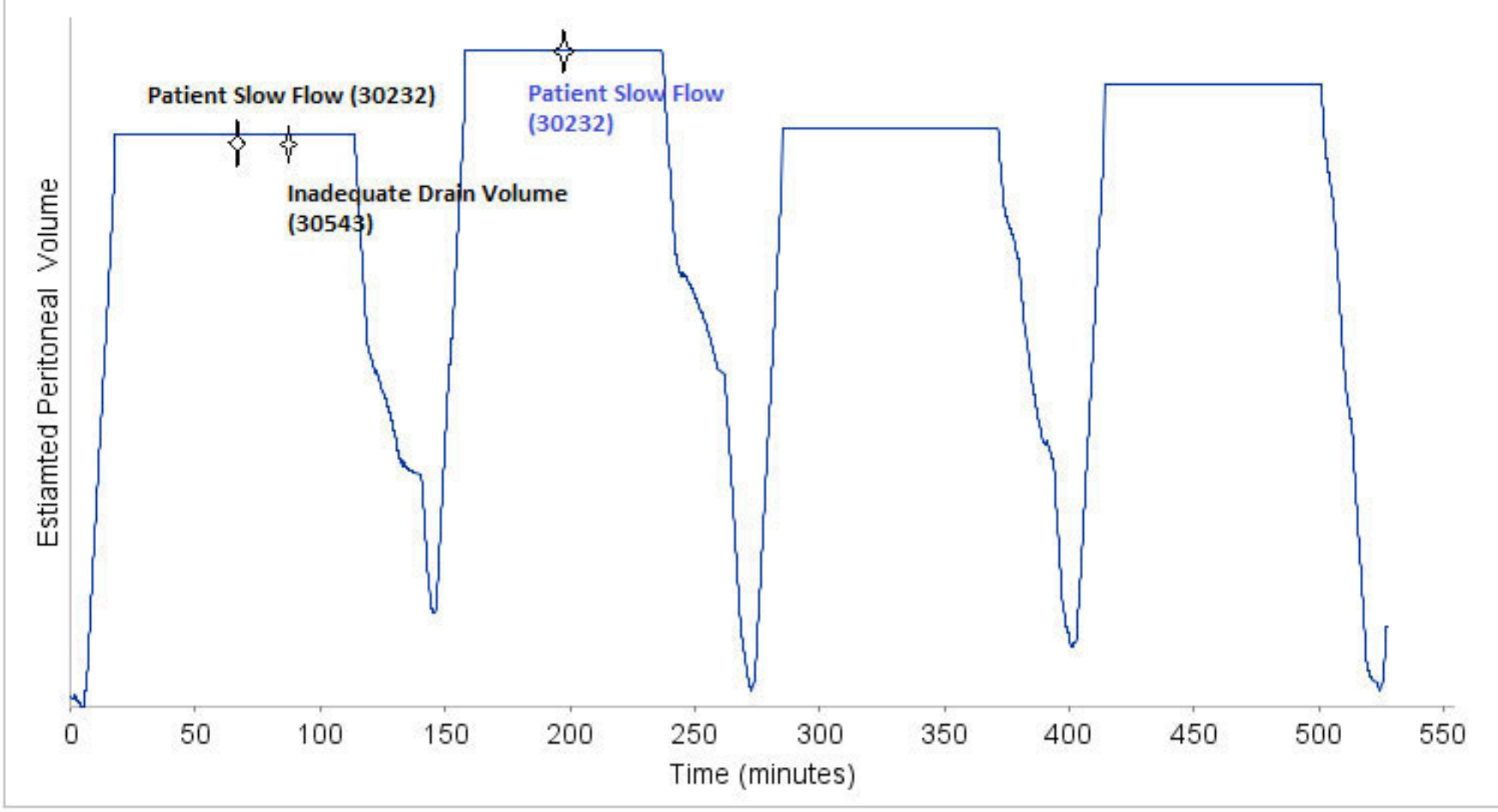

Manual Exchange	Solution Type	Fill Volume	Drain Volume	UF
1	Dianeal	ml	ml	ml
2	Dianeal	ml	ml	ml

Device Program Details

Total Therapy Volume	9,900 mL	Tidal Reserve Volume	--	Extra Last Drain Mode	Always
Last Fill Volume	300 mL	Max Peritoneal Volume	3,600 mL	Extra Drain UF Limit	--
Last Solution Type	Same	Estimated Night UF	1,100 mL	Extra Last Drain Alert	Off
Full Drain Frequency	--	Min Initial Drain Volume	0 mL	Min Initial Drain Time	5 minutes
Tidal Percent	--	Min Night Drain Time	20 minutes	Min Day Drain Volume Percent	--
Tidal Fill Volume	--	Effluent Sample Reminder	Disabled	Min Day Drain Time	--





Device Data, Drug Survey

Device Data		Drugs	Applied
Device ID	27060118	ESA	Not Surveyed
Device Software Version	2.2.0-5	Antibiotics	Not Surveyed
		Heparin	Not Surveyed
		Vitamin D	Not Surveyed

To mark this treatment as reviewed click the Treatment Reviewed button on the right.

Treatment Reviewed

Device Settings

- Patient Actions**
-  Return to Dashboard
-  Claria Snapshot
-  Claria Treatment Summary
-  Claria Device Settings

[<< Previous Patient](#)
[Next Patient >>](#)

Smith, John

Attending Physician: Dr. Barron, Isabell

Date of Birth: 22 November 1977


Clinic Patient ID:

Baxter Patient ID: O-TRAIN-1003042

Homechoice Claria

The settings below apply to the patient's at-home device. Click to View / Edit the corresponding setting, or click on the Add Program link to create a new setting for this patient, up to four programs.

Patient Activation Code: 070-1979-023

 Patient Settings and System Settings cannot be accessed prior to creating a Device Program.

	Name	Modified By	Modified On	Actions
Device Programs				
Add Program	NewTest1	Tyler, Steven	29 March 2018	View Edit Make Active
Add Program	test	Perry, Joe	29 March 2018	View Edit Make Active
Add Program	TestDevice	Whitford, Brad	29 March 2018	View Edit Make Active
Patient Settings				
	Patient Settings	Hamilton, Tom	31 Oct 2017	View Edit
System Settings				
	System Settings	Densmore, John	27 Oct 2017	View Edit

Device Settings
Edit Device Time

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Home Choice Claria

These settings apply to a patient's device. Please enter all of the settings and select Review to continue to the settings confirmation page.

Apply Clinic Template
All

- Time
- Volume
- Tidal
- Settings

Device Program Name
namePlaceholder

Therapy Mode

Night Terapy Time
 Hours Minutes
(Range: 0:10 - 48:00, in increments of 0:10)

Smart Dwells

Next

Cancel

Delete

Device Settings
Edit Device Volume

- Patient Actions**
- Return to Dashboard
 - Claria Snapshot
 - Claria Treatment Summary
 - Claria Device Settings

<< Previous Patient
Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
 Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Home Choice Claria

These settings apply to a patient's device. Please enter all of the settings and select Review to continue to the settings confirmation page.

Apply Clinic Template
 Apply Template

- Time
- Volume
- Tidal
- Settings

Select if treatment includes the following options:

Total Volume (mL)

Range: 200 - 80,000 mL
The Total Volume is the target therapy volume to be delivered to the patient. The Estimated Therapy Volume will be calculated from the fill volumes and tidal settings and may not equal the Total Volume.

Day Therapy

Night Therapy	Night Fill Volume mL	Number of Night Cycles	Night Therapy Volume (mL)
Yes	<input style="width: 80%;" type="text"/>	0	0
	<i>(Range: 100 - 3,000 mL)</i>		

Last Fill

Same Detox

Last Fill Volume mL

(Range: 100 - 3,000 mL)

Next
Previous

Cancel
Delete

Device Settings

Add New Program

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Home Choice Claria

These settings apply to a patient's device. Please enter all of the settings and select Review to continue to the settings confirmation page.

Apply Clinic Template
All

- Time
- Volume
- Tidal
- Settings

Device Program Name

Therapy Mode

Night Terapy Time
 Hours Minutes
(Range: 0:10 - 48:00, in increments of 0:10)

Smart Dwells

Next

Cancel

Device Settings
Program New Device

- Patient Actions
- Return to Dashboard
- Claria Snapshot
- Claria Treatment Summary
- Claria Device Settings

<< Previous Patient Next Patient >>

Smith, John Attending Physician: Dr. Barron, Isabell
Date of Birth: 22 November 1977 Clinic Patient ID: Baxter Patient ID: O-TRAIN-1003042

Home Choice Claria

These settings apply to a patient's device. Please enter all of the settings and select Review to continue to the settings confirmation page.

Apply Clinic Template
All

- Time
- Volume
- Tidal
- Settings

Device Program Name

Therapy Mode

Night Terapy Time
 Hours Minutes
(Range: 0:10 - 48:00, in increments of 0:10)

Smart Dwells

Next

Cancel