

Infusion Dashboard Dev Audit

Audit performed on a Mac, using Chrome at 1366 px wide. [According to research this is the second most common resolution in the world at roughly 14%.](#)

Legend


Bubbles

- 1 Completed Task
- 2 Concern to be addressed
- 3 Rejected Concern
- 4 NPE/Product Owner comments
- 5 Being discussed or researched
- 6 High Priority

Comments

- Completed Task
- Concern to be addressed
- ~~Completed Task~~
- NPE/Product Owner comments
- Being discussed or researched
- (!) High Priority (!)

Nurse View

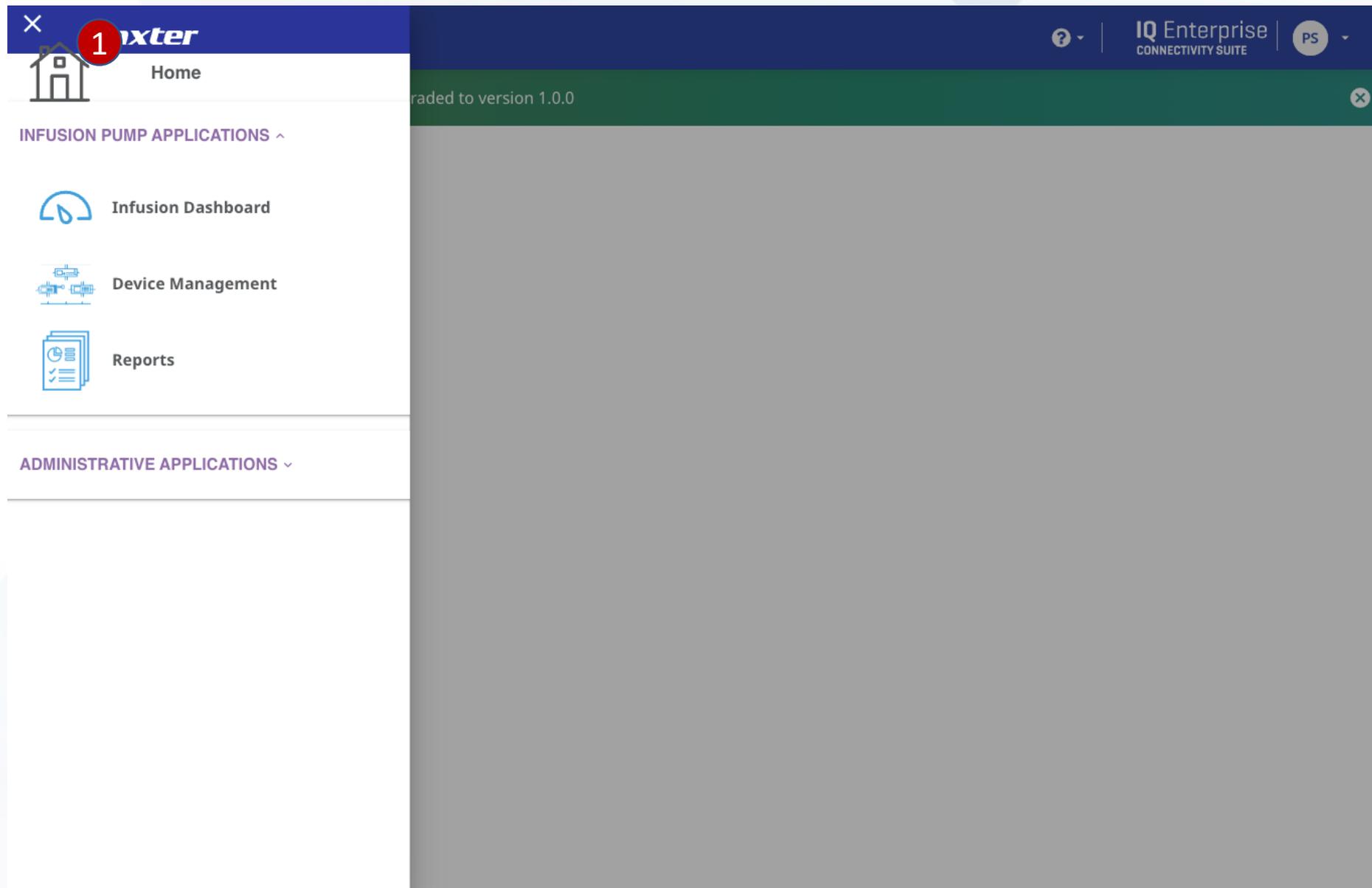

Navigation

1. Home icon above header (from patient index manager)



Home

1

INFUSION PUMP APPLICATIONS ^

- Infusion Dashboard
- Device Management
- Reports

ADMINISTRATIVE APPLICATIONS v

All Filters

Drug Library Name: All | Drug Library Version: All | Device OS: All | Device Status: All

MR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3		Active	19-JUL-2023, 05:36
MR_0265249D48E8-	Novum IQ	SYR		NS-DBG-v1.1.3		Unregistered	20-JUL-2023, 09:13
MR_02BF899A8BD6-	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16

Navigation

1. The blues should align (device management)

× **Baxter**

Infusion Dashboard
SOFTWARE

PS

Home

1

INFUSION PUMP APPLICATIONS ▾

ADMINISTRATIVE APPLICATIONS ▾

SEARCH

[Manage Locations](#)

[Manage Patients](#)

[Manage Devices](#)

[Audit](#)

Infusion Status Legend !

No data available

< >

0 Total items

1.0.10 | [Privacy Policy](#) | [Terms & Conditions](#)

Navigation

1. Fonts aren't available and the overlay covers the top when it doesn't in other areas (dashboard)

Baxter | Infusion Dashboard | PS

Nurse View

Search by Patient Identifier | Enter Patient Identifier | SEARCH | [Manage Locations](#) | [Manage Patients](#) | [Manage Devices](#) | [Audit](#)

0 Patients Alarming | 4 Patients | Infusion Status Legend

Bob,Johnson	MRN-10192	1 Male	Childrens Hospital,Pediatric,6,2	FLOW SHEET →
John,Smith	MRN-10121	0 Male	Childrens Hospital ,NICU,5,1	FLOW SHEET →
Cindy ,Smith	MRN-20121	1 Male	Childrens Hospital,Pediatric,2,2	FLOW SHEET →
Jane,Wilson	MRN -12100	0 Male	Childrens Hospital,NICU,10,2	FLOW SHEET →

25 Items per page | 4 Total items

© 2023 Baxter | 1.0.10 | [Privacy Policy](#) | [Terms & Conditions](#)

Dashboard

1. No fonts are rendering correctly
2. Gray strip has a 1-2 margin from both edges
3. Why is there so much spacing between the buttons? Why are the buttons in the search/filter row? They feel like filters not buttons
4. Light gray makes the field appear inactive
5. Why is the "4" the same color as the text links?
6. Why does clicking to the right of the arrow launch the flow sheet?
7. The footer is too tall, it should be less than ¼ the height

Nurse View

Search by Patient Identifier

Enter Patient Identifier

SEARCH

[Manage Locations](#)

[Manage Patients](#)

[Manage Devices](#)

[Audit](#)

0 Patients Alarming

4 Patients

Infusion Status Legend

Bob,Johnson	MRN-10192	1 Male	Childrens Hospital,Pediatric,6,2	FLOW SHEET →
Insulin(DKA) 50 Units / 100 mL Status Infusing (1) LAST UPDATED 1831 hr(s) 53 min(s) ! (5) TIME REMAINING UNKNOWN (2) VOLUME REMAINING 24.11 mL DOSE RATE 0.4 Units/kg/hr + (6) FLOW RATE 56 mL/hr (3) VOL INFUSED 50.89 mL BAG VOL 75 mL (4)				
John,Smith	MRN-10121	0 Male	Childrens Hospital ,NICU,5,1	FLOW SHEET →
Cindy ,Smith	MRN-20121	1 Male	Childrens Hospital,Pediatric,2,2	FLOW SHEET →
Jane,Wilson	MRN -12100	0 Male	Childrens Hospital,NICU,10,2	FLOW SHEET →

25 Items per page

< 1 >

4 Total items

Dashboard

1. Infusing is not vertically aligned with the play button
2. The horizontal spacing is inconsistent and crowded
3. The bag icon is floating randomly and gets lost.
4. Is the table scroll necessary?
5. Not a dev issue: unless this is a button this needs to be done differently
6. The yellow is difficult to read and the icon should be "inline"

Nurse View

[Manage Locations](#)
[Manage Patients](#)
[Manage Devices](#)
[Audit](#)

Patient Name 1
Patient Identifier 4 Patients
 Infusion Status Legend !

▼	Bob,Johnson	MRN-10192	1 Male	Childrens Hospital,Pediatric,6,2	FLOW SHEET →
Insulin(DKA) 50 Units / 100 mL					
<small>Status</small> ▶ Infusing					
<small>LAST UPDATED</small> 1832 hr(s) 01 min(s) !					
<small>TIME REMAINING</small> UNKNOWN					
<small>VOLUME REMAINING</small> 24.11 mL					
<small>DOSE RATE</small> 0.4 Units/kg/hr					
<small>FLOW RATE</small> 56 mL/hr					
<small>VOL INFUSED</small> 50.89 mL					
<small>BAG VOL</small> 75 mL					
^	John,Smith	MRN-10121	0 Male	Childrens Hospital ,NICU,5,1	FLOW SHEET →
^	Cindy ,Smith	MRN-20121	1 Male	Childrens Hospital,Pediatric,2,2	FLOW SHEET →
^	Jane,Wilson	MRN -12100	0 Male	Childrens Hospital,NICU,10,2	FLOW SHEET →

Dashboard

1. Dropdown should not cover form field. Why is "Patient Identifier" pre-selected?

Nurse View

Search by Patient Identifier

Enter Patient Identifier

SEARCH

[Manage Locations](#)

0 Patients Alarming

4 Patients

Bob,Johnson

MRN-10192

1 Male

Childrens Hosp

Insulin(DKA) 50 Units / 100 mL

Status

Infusing

LAST UPDATED

1833 hr(s) 41 min(s) !

TIME REMAINING

UNKNOWN

VOLUME REMAINING

24.11 mL

John,Smith

MRN-10121

0 Male

Childrens Hosp

Cindy ,Smith

MRN-20121

1 Male

Childrens Hosp

Jane,Wilson

MRN -12100

0 Male

Childrens Hosp

25

Items per page

<

1

>

Infusion Status

The Baxter Infusion Dashboard provides an at-a-glance view of assigned patients and the status of their infusions.

Patients experiencing an alarm will rise to the top of the list. Below the alarm icon, an additional stop or play icon will indicate whether or not the infusion is still running.



High Priority Alarm: red, triangle with three exclamation points. Infusions will usually automatically STOP, but not always



Medium and Low Priority Alarms: low (one exclamation point yellow), medium (two exclamation points yellow). Infusions will continue running



Infusing: green, play icon. Indicates infusion is running as expected..



Paused: gray, pause icon. Indicates infusion has been temporarily stopped and can be resumed.



Completed: blue, check icon. Indicates an infusion that has concluded and been closed by a clinician in the past 24 hours.

The following icons indicate the type of device associated with the infusion



Large Volume Pump (LV)



Syringe Pump

Dashboard

1. Missing X to close overlay
2. If you click the background where a Flow Sheet button is located it navigates to flow sheet
3. Margins are too tight
4. Vertical spacing doesn't make it obvious the title refers to the icons. Icon labels shouldn't be bolded compared to the title.

Manage Locations

1 Clinician: Peter S

Available Locations

My Locations

Search by Location or Facility or Department

Search by Location or Facility or Department

Location ID	Facility	Department
<input type="checkbox"/>	Childrens Hospital ,NICU	Childrens Hospital
<input type="checkbox"/>	Childrens Hospital,Emergency	Childrens Hospital
<input type="checkbox"/>	Childrens Hospital,NICU	Childrens Hospital
<input type="checkbox"/>	Childrens Hospital,Oncology	Childrens Hospital
<input type="checkbox"/>	Childrens Hospital,Pediatric	Childrens Hospital

Location ID	Facility	Department
No data available		

3

ADD

REMOVE

SAVE CHANGES

CANCEL

Dashboard: Manage Locations

1. Physician name is floating
2. If user can only select one row this should be radio buttons not check boxes
3. Why are the buttons blue when "off" but white with a blue outline when "on?"
4. The spacing on the commas is not aligned in the rows

Manage Locations

Clinician: Peter S



Available Locations

Search by Location or Facility or Department

My Locations

Search by Location or Facility or Department

Location ID Facility Department

1	Childrens Hospital,Emergency	Childrens Hospital
	Childrens Hospital,NICU	Childrens Hospital
	Childrens Hospital,Oncology	Childrens Hospital
	Childrens Hospital,Pediatric	Childrens Hospital

Location ID Facility Department

<input type="checkbox"/>	Childrens Hospital ,NICU	Childrens Hospital
--------------------------	--------------------------	--------------------

ADD

REMOVE

SAVE CHANGES

CANCEL

Dashboard: Manage Locations

1. Why can't any other locations be added after one?
2. Does Facility have a sort or is it not available because the facilities are the same?

Manage Locations

Clinician: Peter S

Available Locations

My Locations

Search by Location or Facility or Department

Search by Location or Facility or Department

Location ID	Facility	Department	Location ID	Facility	Department
Childrens Hospital,Pediatric	Childrens Hospital				
Childrens Hospital,Oncology	Childrens Hospital				
Childrens Hospital,NICU	Childrens Hospital				
Childrens Hospital,Emergency	Childrens Hospital				

Confirm 1

Are you sure you want to save this changes ?

3 CANCEL 2 OK

ADD

REMOVE

SAVE CHANGES

CANCEL

Dashboard: Manage Locations

1. What is the top label supposed to look like? The font is too light
2. The button should be a blue pill
3. Cancel is too light

Nurse View



Search by Patient Identifier

Enter Patient Identifier

SEARCH

[Manage Locations](#)

[Manage Patients](#)

[Manage Devices](#)

[Audit](#)

0 Patients Alarming

4 Patients

Infusion Status Legend

Bob,Johnson MRN-10192 1 Male Childrens Hospital,Pediatric,6,2 FLOW SHEET →

Insulin(DKA) 50 Units / 100 mL

Status
▶ Infusing

LAST UPDATED
1833 hr(s) 41 min(s) !

TIME REMAINING
UNKNOWN

VOLUME REMAINING
24.11 mL

DOSE RATE
0.4 Units/kg/hr
⬆

FLOW RATE
56 mL/hr

VOL INFUSED
50.89 mL

BAG VOL
75 mL

John,Smith MRN-10121 0 Male Childrens Hospital ,NICU,5,1 FLOW SHEET →

DOBUtamine 500 mg / 250 mL

Status
▶ Infusing

LAST UPDATED
1833 hr(s) 40 min(s) !

TIME REMAINING
UNKNOWN

VOLUME REMAINING
41 mL

DOSE RATE
40 mcg/kg/min
⬆

FLOW RATE
84 mL/hr

VOL INFUSED
71.65 mL

BAG VOL
100 mL

Cindy,Smith MRN-20121 1 Male Childrens Hospital,Pediatric,2,2 FLOW SHEET →

25 Items per page

< 1 >

4 Total items

Dashboard: Manage Locations Added

1. The badge has opacity and should be 100%, also this is different than other badges

Nurse View

Search by Patient Identifier 

Enter Patient Identifier

SEARCH

[Manage Locations](#)

[Manage Patients](#)

[Manage Devices](#)

[Audit](#)

0 Patients Alarming 

1 Patients

Infusion Status Legend 

	Jane, Wilson	MRN - 12100	0 Male	Childrens Hospital, NICU, 10, 2	FLOW SHEET 
---	--------------	-------------	--------	---------------------------------	--

25 

Items per page

1 Total items

Dashboard

1. After making any patient changes I lost all the patients but one

Manage Patients

Clinician: Peter S

1



Available Patients

Patient Identifier ▾

Enter Patient Identifier or ...

SEARCH

2

My Patients

Enter Patient Identifier or Patient Name

0 Patients

<input type="checkbox"/>	Patient Identifier ▾	Patient Name	Age/Gender	Location ID
<input type="checkbox"/>	MRN -12100	Jane, Wilson	0 / M	Childrens Hospital, NICU, 10, 2

4

3

<input type="checkbox"/>	Patient Identifier ▾	Patient Name	Age/Gender	Location ID
No data available				

ADD ▶

◀ REMOVE

SAVE CHANGES

CANCEL

Dashboard: Manage Patients

1. Dr name is "floating"
2. Search button is not lined up to form field
3. Commas are off
4. Could not click in row checkbox

< BACK TO DASHBOARD

1

Patient Device Association

Search by Patient Identifier

Patient Identifier	Patient Name	Age/Gender	Facility / Department	Room / Bed	No. Of Associated Devices	Manage Devices
MRN -12100	Jane,Wilson	0 / M	Childrens Hospital /	10 / 2	2	

10 Items per page 1 Total items

Dashboard: Patient Device Association

- 1. Fonts are not correct

[← BACK TO DASHBOARD](#)

Patient Device Association

Patient Identifier

1

Patient Name

Age/Gender

Facility / Department

Room / Bed

No. Of Associated Devices

Manage Devices

MRN -12100

Jane, Wilson

0 / M

Childrens Hospital /

10 / 2

2



10

Items per page



1



1 Total items

Dashboard: Patient Device Association

1. Dropdown shouldn't cover field

BACK TO DASHBOARD

Audit

1

2023-10-03

3

EXPORT AS CSV

Context	Event	User Information	Timestamp	Old Values	New Values
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:26	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN-10121	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:26	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN-20121	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:26	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN -12100	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:26	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN-10192	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:26	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN-2913	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:27:48	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN -12100	
Nurse to Patient	Deleted	Name: nurse User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 IP Address: 192.168.254.1	03 Oct 2023 16:39:06	User Id: 75ec0dd9-e639-4d6f-9c95-4b1f0f2a00a7 Patient Identifier: MRN -12100	

4

10

Items per page

< 1 >

9 Total items

Dashboard: Audit

- Calendar icon does not work
- User Information, Time Stamp and Old Values do not sort
- Should the button be at the bottom?
- User Information is difficult to scan

Device Management

1

Device List 6

Advanced Filters

<input type="checkbox"/> 2	Serial Number	Hospital Device ID 3	EMR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
<input type="checkbox"/>	0265249D48E8-1		deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3		Active	19-JUL-2023, 05:36
<input type="checkbox"/>	0265249D48E8-5		EMR_0265249D48E8-5	Novum IQ	SYR		NS-DBG-v1.1.3		Unregistered	20-JUL-2023, 09:13
<input type="checkbox"/>	02BF899A8BD6-1 4		EMR_02BF899A8BD6-1	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-2		deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-3		deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-4		deviceEmrId	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-5		deviceEmrId5	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	06DCE2723FAC-46		deviceEmrId46	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-47		deviceEmrId47	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-48		deviceEmrId48	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38

10

Items per page 5

1

< 1 2 >

12 Total items

EXPORT LIST

Device Management

1. Navigating to screen the header and footer do not appear
2. What does triple dot overflow menu do?
3. Column names are not inline
4. "Word wrap" should not happen, it is better to have a horizontal scroll. This makes scanning more difficult
5. Items per page should be inline
6. Why is there a tab if there are no other options to select?

Device Management

Device List

Advanced Filters ▾

<input type="checkbox"/>	Serial Number	Hospital Device ID	EMR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
<input type="checkbox"/>	06DCE2723FAC-50		deviceEmrId50	Novum IQ	SYR				Unregistered	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-49		deviceEmrId49	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-48		deviceEmrId48	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-47		deviceEmrId47	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	06DCE2723FAC-46		deviceEmrId46	Novum IQ	SYR	NPP			Active	19-JUL-2023, 05:38
<input type="checkbox"/>	02BF899A8BD6-5		deviceEmrId5	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-4		deviceEmrId	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-3		deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-2		deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-1		EMR_02BF899A8BD6-1	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16

10 ▾ Items per page < 1 2 > 12 Total items EXPORT LIST

Device Management

- 1. Sort arrow should be inline with the column name, not below

Device Management

Device List

Advanced Filters ^

[Reset All Filters](#) 2

Search By 5

Device Model Device Type Drug Library Name Drug Library Version Device OS Device Status

All All 4 All All All All

APPLY FILTERS 3 ENTERPRISE GROUP

<input type="checkbox"/>	:	Serial Number	Hospital Device ID	EMR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
<input type="checkbox"/>		0265249D48E8-1		deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3		Active	19-JUL-2023, 05:36
<input type="checkbox"/>		0265249D48E8-5		EMR_0265249D48E8-5	Novum IQ	SYR		NS-DBG-v1.1.3		Unregistered	20-JUL-2023, 09:13
<input type="checkbox"/>		02BF899A8BD6-1		EMR_02BF899A8BD6-1	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-2		deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-3		deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16

10 Items per page 1

< 1 2 >

12 Total items

EXPORT LIST

Device Management

1. Top and bottom appeared when Advanced Filters was clicked
2. Reset Filters is in a "weird" place
3. Apply Filters should be disabled until filters have been selected
4. It is not clear all filter selections are disabled until a search is selected
5. How does the search work? Where is the search button?

Device Management

Device List

Advanced Filters ^

[Reset All Filters](#)

Serial Number **1**

Hospital Device ID

EMR ID

APPLY FILTERS ENTERPRISE GROUP

Drug Library Name: All | Drug Library Version: All | Device OS: All | Device Status: All

<input type="checkbox"/>	:	Serial Number	Hospital Device ID	EMR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
<input type="checkbox"/>		0265249D48E8-1		deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3		Active	19-JUL-2023, 05:36
<input type="checkbox"/>		0265249D48E8-5		EMR_0265249D48E8-5	Novum IQ	SYR		NS-DBG-v1.1.3		Unregistered	20-JUL-2023, 09:13
<input type="checkbox"/>		02BF899A8BD6-1		EMR_02BF899A8BD6-1	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-2		deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-3		deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16

Device Management

1. Dropdown should not be over form field. Why is Hospital Device ID available when there is no data in the column?

Device Management

Device List

Advanced Filters ^

[Reset All Filters](#)

Search By EMR ID

Device Model: All | Device Type: All | Drug Library Name: All | Drug Library Version: All | Device OS: All | Device Status: All

APPLY FILTERS 1 **ENTERPRISE GROUP**

EMR ID: device x 2

<input type="checkbox"/>	:	Serial Number	Hospital Device ID	EMR ID	Device Model	Device Type	Enterprise Group	Device Firmware	Library Name (Version)	Device Status	Last Connection
<input type="checkbox"/>		0265249D48E8-1		deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3		Active	19-JUL-2023, 05:36
<input type="checkbox"/>		02BF899A8BD6-2		deviceEmrId2	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-3		deviceEmrId3	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-4		deviceEmrId	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16
<input type="checkbox"/>		02BF899A8BD6-5		deviceEmrId5	Novum IQ	LVP	NPP			Active	19-JUL-2023, 06:16

Items

Device Management

1. Why does Apply Filters function as a search button? Why don't any of the filter dropdowns work?
2. Why can only one filter be applied at a time? I tried to do a device and serial number as a combined filter and it only works for the latest filter.

Device Management

Device List

Advanced Filters ^

[Reset All Filters](#)

Search By

Device Model	Device Type	Drug Library Name	Drug Library Version	Device OS	Device Status
All					

APPLY FILTERS

Enterprise Hierarchy

Please Select Enterprise Group **1**

APPLY **2** [Cancel](#)

<input type="checkbox"/>	Serial N	Device Model	Device Type	Drug Library Name	Drug Library Version	Device OS	Device Status	Last Connection
<input type="checkbox"/>	0265249D48E8-1	deviceEmrId1	Novum IQ	LVP	NPP	NS-DBG-v1.1.3	Active	19-JUL-2023, 05:36
<input type="checkbox"/>	0265249D48E8-5	EMR_0265249D48E8-5	Novum IQ	SYR		NS-DBG-v1.1.3	Unregistered	20-JUL-2023, 09:13
<input type="checkbox"/>	02BF899A8BD6-1	EMR_02BF899A8BD6-1	Novum IQ	LVP	NPP		Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-2	deviceEmrId2	Novum IQ	LVP	NPP		Active	19-JUL-2023, 06:16
<input type="checkbox"/>	02BF899A8BD6-3	deviceEmrId3	Novum IQ	LVP	NPP		Active	19-JUL-2023, 06:16

Device Management

1. How do I access the enterprise groups? If one isn't available why isn't there a message saying that?
2. Too much space between Cancel and Apply.

Reports

1

Library [Select a Drug Library to View Reports](#)

Quick Search

Advanced Filters

- Compliance**
- Limits
- Dose Rate Change
- Device Usage

[Select a Drug Library to View Reports](#)

1

Reports

1. Header and footer missing when first navigating to page
2. Why are dates available before selecting drug reports?
3. Apply button is too far away from everything
4. Why are advanced filters unavailable?

Reports

Select a Drug Library to View Reports

Select Drug Library ✕

Search By Library Name SEARCH

Library Name	Device Model	Device Type	Imported By	Import Date
<input type="checkbox"/> 1 novumLVPSYRUS_Demo	NOVUM IQ	LVP	ap	19-JUL-2023
<input type="checkbox"/> npp	NOVUM IQ	LVP	ap	19-JUL-2023

10 Items per page **2** < 1 > 2 items VIEW REPORTS

Reports

- 1. Since only one report can be selected it should be radio buttons not checkmarks
- 2. Should be inline

Reports

Select Drug Library ×

1 Search By Library Name × SEARCH

2 Library Name: novum ×

Library Name	Device Model	Device Type	Imported By	Import Date
<input type="checkbox"/> NovumLVPSYRUS_Demo	NOVUM IQ	LVP	ap	19-JUL-2023

10 ▾ Items per page < 1 > 1 items VIEW REPORTS

Reports

- 1. Why does the field say “Search by Library Name?” when the device model is the search result?
- 2. Why is there a “filter” bubble when this is a search?

Reports

1. Change drug library is disconnected from drug library
2. The dates in blue look like buttons but are not clickable. When I change the dates in the text link it is not reflected in the blue (note: text link shows oct 3 while blue shows oct 2).
3. What is a quick search? It functions like a filter.
4. Clear filters is separated from everything. Also should be below advanced filters as I am assuming it will clear all filters. (it's not obvious when filters are cleared as nothing changes)
5. See above (4)
6. It is not clear checks control doughnut chart. It is not clear they are related
7. Gray is difficult to read. Also, are there other colors in that area?
8. When exporting pdf why is there an export button in the pdf?

Reports

Library: NovumLVPSYRUS_Demo [1 Change Drug Library](#)

Quick Search [3](#) Search By [01-JAN-2023 - 03-OCT-2023](#) [2](#) [APPLY](#)

Start Date: [01-JAN-2023](#) End Date: [02-OCT-2023](#) [2](#)

Advanced Filters [4 Clear Filters](#) [5 Save / Manage Filters](#)

Compliance Limits Dose Rate Change Device Usage

View Infusions:

<input checked="" type="checkbox"/>	Auto Program	0	0.00%
<input checked="" type="checkbox"/>	Manual	1	100.00%
<input checked="" type="checkbox"/>	DERS	1	100.00%
<input checked="" type="checkbox"/>	Basic Mode	0	0.00%
<input checked="" type="checkbox"/>	Manual + Basic	0	0.00%
<input checked="" type="checkbox"/>	Auto + Basic	0	0.00%
<input checked="" type="checkbox"/>	Auto + DERS	0	0.00%
<input checked="" type="checkbox"/>	Manual + DERS	1	100.00%

[6](#) [EXPORT](#) [8](#)

[7](#) [EXPORT](#)

Shifts **1** Shift 1 - 07:00 to 15:00 Shift 2 - 15:00 to 23:00 Shift 3 - 23:00 to 07:00

Enterprise Units Default NPP **3**

- Care Areas **2**
- Adult ICU
 - Adult ICU AP
 - Anesthesia
 - Emergency
 - L&D
 - MedSurg
 - NICU
 - NICU AP

Drugs **2**

Enter a name

APPLY

Compliance Limits Dose Rate Change Device Usage

View Infusions:			
<input checked="" type="checkbox"/>	Auto Program	0	0.00%
	Manual	1	100.00%




Reports

- Shifts, Enterprise Units, Care Area & Drugs should at least be bold or a larger font to make more clear these are labels to form a hierarchy
- Are Care Areas and Drugs related?
- Are there other possible options? Is this an on/off? If yes it should be a toggle.

View Infusions:

<input checked="" type="checkbox"/>	Auto Program	0	0.00%
	Manual	1	100.00%
<hr/>			
<input checked="" type="checkbox"/>	DERS	1	100.00%
	Basic Mode	0	0.00%
<hr/>			
<input checked="" type="checkbox"/>	Manual + Basic	0	0.00%
	Auto + Basic	0	0.00%
	Auto + DERS	0	0.00%
	Manual + DERS	1	100.00%

EXPORT

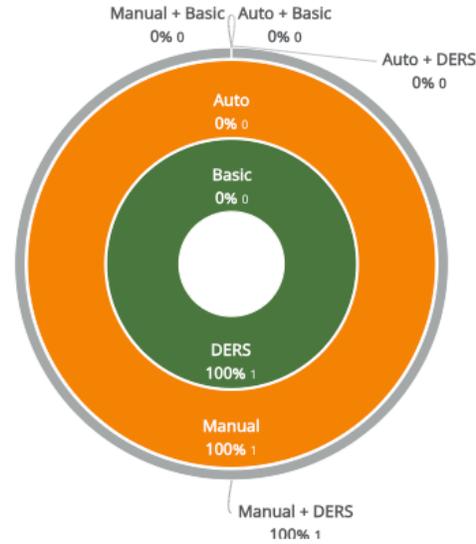

EXPORT

Infusions by Care Area

2

Care Area	Total Infusions	Auto Program	Manual	DERS	Basic Mode	Manual + Basic	Auto + Basic	Auto + DERS	Manual + DERS
All Selected Care Areas	1	0 0.00%	1 100.00%	1 100.00%	0 0.00%	0 0.00%	0 0.00%	0 0.00%	1 100.00%
Adult ICU	1	0 0.00%	1 100.00%	1 100.00%	0 0.00%	0 0.00%	0 0.00%	0 0.00%	1 100.00%

EXPORT

Reports

1. Header re-appeared when I logged back in
2. When clicking Infusion by Care Area it changed the sorting of All Selected Care Area
3. There should be no word wrap in columns, horizontal scroll is better
4. The text links are confusing

Reports

Library

Quick Search

Start Date: 01-JAN-2023

Advanced Filters

Compliance Limits Dose Rate

View Infusions:

- Auto Program
Manual
- DERS
Basic Mode
- Manual + Basic
Auto + Basic
Auto + DERS
Manual + DERS

Select a Date Range

Time Frame: Custom

From: 2023 Sun, Jan 1

January 2023						
S	M	T	W	T	F	S
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

To: 2023 Mon, Oct 2

October 2023						
S	M	T	W	T	F	S
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

OK **1** Cancel

EXPORT

Manual + DERS
100% 1

EXPORT

Reports

- Buttons too far apart

Reports

Library

Quick Search

Start Date: 01-JAN-2023

End Date

Advanced Filters

Compliance Limits Dose Rate

View Infusions:

- Auto Program Manual
- DERS Basic Mode
- Manual + Basic
Auto + Basic
Auto + DERS
Manual + DERS

Select a Date Range

Time Frame

From 2023 **Sun, Jan 1** To 2023 **Mon, Oct 2**

Year to Date
Custom

January 2023						
S	M	T	W	T	F	S
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

October 2023						
S	M	T	W	T	F	S
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30	31				

OK Cancel

[Change Drug Library](#)

APPLY

[Clear Filters](#)

[Save / Manage Filters](#)

EXPORT

Manual + DERS
100% 1

EXPORT

Reports

1. Drop down should not cover form field

Reports

Library

NovumLVPSYRUS_Demo

[Change Drug Library](#)

Quick Search

Search By

[01-JAN-2023 - 03-OCT-2023](#)

APPLY

Start Date: 01-JAN-2023

End Date: 03-OCT-2023

[Clear Filters](#)

Advanced Filters

[Save / Manage Filters](#)

Compliance **2**

Units

Dose Rate Change

Device Usage

[< BACK TO COMPLIANCE REPORT](#)

All Selected Care Areas **1**

3

Date / Time	Drug Name (Modifier)	Concentration	Program Type	Mode	Care Area	Starting Dose Rate
-------------	----------------------	---------------	--------------	------	-----------	--------------------

▼

19-JUL-2023, 05:03:21

ceFAZolin

1 g / 50 mL

Manual

DERS

Adult ICU

90.909 mL/hr

10

Items Per Page



1



1 Items

EXPORT

Reports

1. Looks like a button yet isn't clickable
2. Why is the doughnut gone when clicking a button in a table below the doughnut? That is confusing
3. Text link is wrong color

Reports

Library: NovumLVPSYRUS_Demo [Change Drug Library](#)

Quick Search: Search By [01-JAN-2023 - 02-OCT-2023](#) [APPLY](#)

Start Date: **01-JAN-2023** End Date: **03-OCT-2023** [Clear Filters](#)

Advanced Filters [Save / Manage Filters](#)

Compliance Limits Dose Rate Change Device Usage

[< BACK TO COMPLIANCE REPORT](#) [All Selected Care Areas](#)

Date / Time	Drug Name (Modifier)	Concentration	Program Type	Mode	Care Area	Starting Dose Rate
19-JUL-2023, 00:03:21	ceFAZolin	1 g / 50 mL	Manual	DERS	Adult ICU	90.909 mL/hr
HOSPITAL DEVICE ID	SERIAL NUMBER	LOCATION / ACCESS POINT	DEVICE MODEL & TYPE	INFUSION ID	EMR DEVICE ID	ENTERPRISE GROUP
	02BF899A8BD6-1		Novum IQ LVP	000000001689747345297204600000	EMR_02BF899A8BD6-1	NPP

10 Items Per Page [<](#) [1](#) [>](#) 1 Items [EXPORT](#)

Reports


1. When table row expanded columns shift right
2. Expanded row color should remain the same color as its master row
3. Data rows should be inline

Show: 1

View: Auto Programming Initial Programming Soft Limits
 Manual 2 During Infusion Hard Limits

QUANTITY PERCENTAGE 3

Limit Events Leaderboard

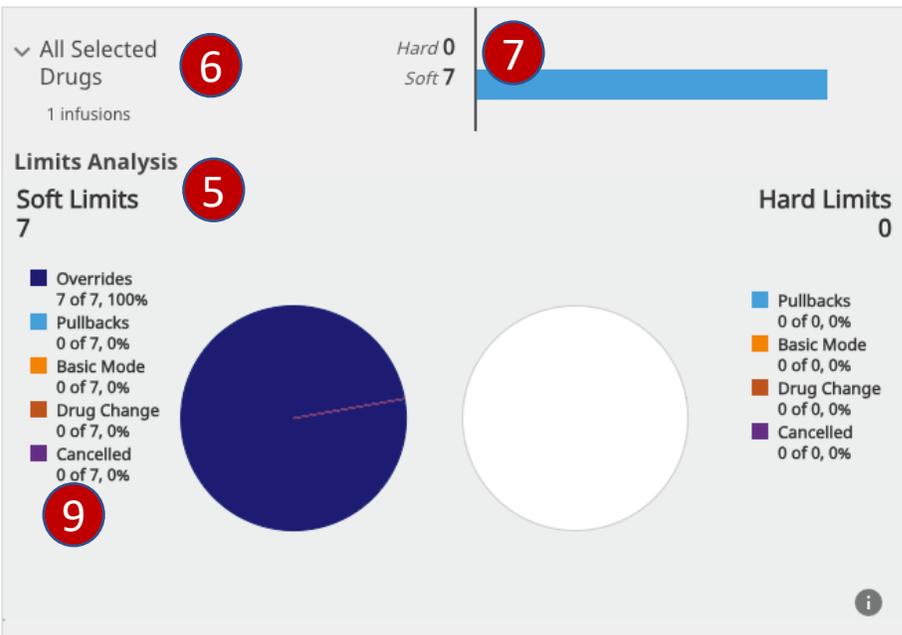

Care Area	Limit Type	Event
<input checked="" type="checkbox"/> 4 All Selected Care Areas	All Limits	7
<input type="checkbox"/> Adult ICU	All Limits	7



Reports

1. The dropdown should not be so tall
2. It's not clear these check boxes control the events
3. It's not clear this toggle controls the views on the charts
4. These should be radio buttons not check marks since user can only select one. Not clear these control the pie chart.
5. Limits Analysis, Soft Limits and Hard Limits: no clear hierarchy due to font sizes
6. Should be inline
7. Should this be below "All Drugs Selected?"
8. Is the bar chart part of the expansion? Is it at the same level or a subset?
9. It's not clear these items control what is visible in the pie chart. I discovered it by accident.

reports

Library: NovumLVPSYRUS_Demo [Change Drug Library](#)

Quick Search: Search By [01-JAN-2023 - 02-OCT-2023](#) [APPLY](#)

Start Date: [01-JAN-2023](#) End Date: [02-OCT-2023](#) [Clear Filters](#)


Advanced Filters [Save / Manage Filters](#)

- Compliance
- Limits**
- Dose Rate Change
- Device Usage

Show: View: Auto Programming Initial Programming Soft Limits
 Manual During Infusion Hard Limits

- QUANTITY**
- PERCENTAGE

Limit Events Leaderboard

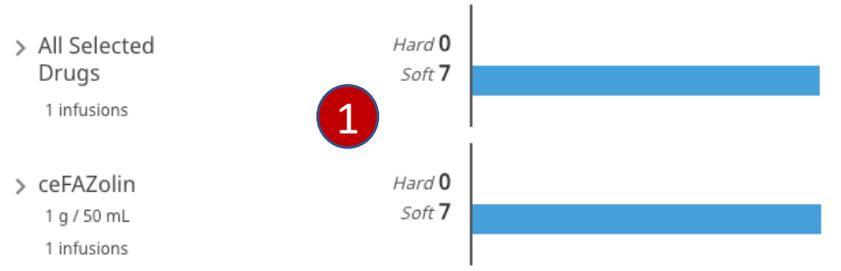

Reports

1. Not enough space between groups. No clear divisions

Start Date: 01-JAN-2023

End Date: 02-OCT-2023

[Clear Filters](#)

Advanced Filters ▾

[Save / Manage Filters](#)

Reports

1. Filters look like buttons
2. Wrong color for text link
3. Columns not inline
4. All rows should be inline

Compliance **Limits** Dose Rate Change Device Usage

Show:

All Limits ▾

View:

Auto Programming

Initial Programming

Soft Limits

Manual

During Infusion

Hard Limits

2 [GO TO LIMITS EVENTS LEADERBOARD](#)

[All Selected Drugs](#) **1** [All Selected Care Areas](#)

Date / Time	Drug Name (Modifier)	Concentration	Limit Type	Alert 3	Next Action	Care Area	Attempted Value	Limits Exceeded
19-JUL-2023, 00:50:17	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	56 min	Soft: 10 - 30 min
19-JUL-2023, 00:50:00 4	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	80 min	Soft: 10 - 30 min
19-JUL-2023, 00:28:14	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	57 min	Soft: 10 - 30 min
19-JUL-2023, 00:24:58	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	82 min	Soft: 10 - 30 min
19-JUL-2023, 00:11:01	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	64 min	Soft: 10 - 30 min
19-JUL-2023, 00:09:07	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	31 min	Soft: 10 - 30 min
19-JUL-2023, 00:03:32	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	33 min	Soft: 10 - 30 min

10 ▾

Items Per Page

< **1** >

7 items

EXPORT

Show: All Limits View: Auto Programming Initial Programming Soft Limits
 Manual During Infusion Hard Limits

[< BACK TO LIMITS EVENTS LEADERBOARD](#)

[All Selected Drugs](#)

[All Selected Care Areas](#)

Date / Time	Drug Name (Modifier)	Concentration	Limit Type	Alert	Next Action	Care Area	Attempted Value	Limits Exceeded
19-JUL-2023, 00:50:17	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	56 min	Soft: 10 - 30 min
<div style="display: flex; justify-content: space-between;"> <div style="width: 30%;"> <p>HOSPITAL DEVICE ID</p> <p>DEVICE SERIAL NUMBER 1</p> <p>02BF899A8BD6-1</p> <p>LOCATION / ACCESS POINT</p> <p>Novum IQ LVP</p> </div> <div style="width: 30%;"> <p>INFUSION ID</p> <p>000000001689747345297204600000</p> <p>DEVICE MODEL & TYPE 2</p> <p>Novum IQ LVP</p> <p>ENTERPRISE GROUP</p> <p>NPP</p> </div> <div style="width: 30%; border: 1px solid #ccc; padding: 5px;"> <p>LIMIT EVENT PARAMETERS View Infusion Story 4</p> <p>EXISTING DOSE RATE 3 EXISTING VALUE</p> <p>22 mL/hr</p> </div> </div>								
19-JUL-2023, 00:50:00	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	80 min	Soft: 10 - 30 min
19-JUL-2023, 00:28:14	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	57 min	Soft: 10 - 30 min
19-JUL-2023, 00:24:58	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	82 min	Soft: 10 - 30 min
19-JUL-2023, 00:11:01	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	64 min	Soft: 10 - 30 min
19-JUL-2023, 00:09:07	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	31 min	Soft: 10 - 30 min
19-JUL-2023, 00:03:32	ceFAZolin	1 g / 50 mL	Time	DERS - Soft	Override	Adult ICU	33 min	Soft: 10 - 30 min

10 Items Per Page

< 1 >

7 items

[EXPORT](#)

Reports

- Expanded row background should be same color throughout
- Too much vertical space between headers and data
- Very difficult to distinguish between sub square and rest of row due to subtle color variations
- Should be inline

INFUSION STORY SOFTWARE

1



< PREVIOUS INFUSION

DEVICE SERIAL NUMBER 02BF899A8BD6-1

NEXT INFUSION >

INFUSION #

000000001689747345297204600000

INFUSION START DATE & TIME

19-JUL-2023, 00:03:21

INFUSION END DATE & TIME

19-JUL-2023, 01:04:08

CARE AREA

Adult ICU

2

DRUG NAME (MODIFIER)

ceFAZolin

CONCENTRATION

1 g / 50 mL

STARTING DOSE RATE

90.909 mL/hr

STARTING FLOW RATE

90.909 mL/hr

DEVICE MODEL & TYPE

Novum IQ LVP

ENTERPRISE GROUP

NPP

DRUG ALIAS

cefazolin


PROGRAM TYPE

Manual

3

TOTAL INFUSION TIME 1 hour 47 seconds

TOTAL VOLUME INFUSED 50 mL

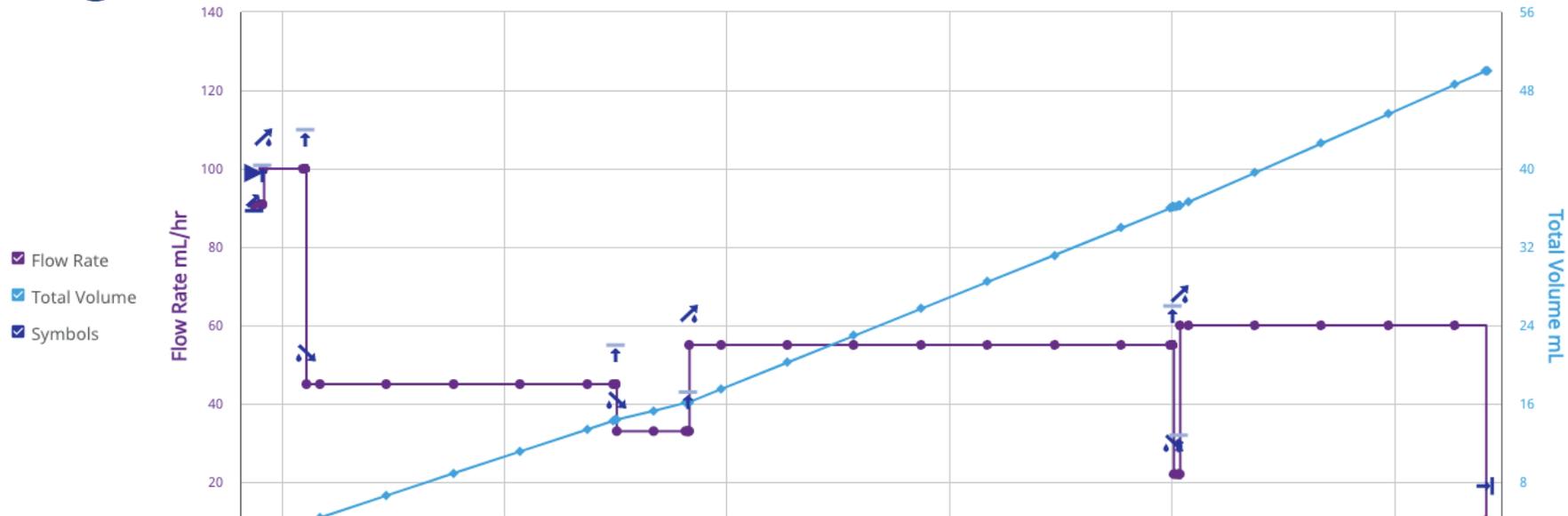

Reports

1. No clear hierarchy with font sizes
2. Too much vertical spacing between label and data. Too little vertical spacing between rows
3. Icon is floating and has no context

Change Type: **1** All Dose Rate Changes

Sort By: **1** Total Number of Dose Rate Changes

2 All Selected Drugs

100% infusions had at least one dose rate change 1 of 1 infusions

3

6 total dose rate changes over 1 infusions

4 0 (0%) Overrides 6 (100%) Within Limits

5

7 ceFAZolin

1 g / 50 mL

6

Dose Change Programming Mode - All Selected Drugs

<input checked="" type="checkbox"/>	Auto Program	0	0%
<input checked="" type="checkbox"/>	DERS	6	100%
<input checked="" type="checkbox"/>	Basic Mode	0	0%
<input type="checkbox"/>	Manual	0	0%

Care Area	Programming Mode	Dose Change
<input checked="" type="checkbox"/> All Selected Care Areas	All	6
<input type="checkbox"/> Adult ICU	All	6

8

Reports

1. Too much space between label and field. It appears as Change Type controls the left column and Sort By controls the right column
2. No clear hierarchy based on font sizes.
3. Too much vertical space between label and chart
4. Green and red fonts difficult to read
5. Arrow too subtle and it's in an odd place
6. Too much vertical space
7. It is not clear to click title to see additional details
8. If single select should be radio buttons. It's not clear checkmarks control graphs

Change Type: All Dose Rate Changes

Sort By: Total Number of Dose Rate Changes

All Selected Drugs
100% infusions had at least one dose rate change 1 of 1 infusions



6 total dose rate changes over 1 infusions



1

Dose Change Rate: Average and Ranges:

Limits	Pullback- % attempts	Limit Override- % accepted	Within limits- % accepted
0%	0% (0% to 0%)	0% (0% to 0%)	83.1% (10% to 172.7%)
0%	0% (0% to 0%)	0% (0% to 0%)	-47.2% (-60% to -26.7%)


2

ceFAZolin
1 g / 50 mL



Dose Change Programming Mode - All Selected Drugs

<input checked="" type="checkbox"/>	Auto Program	0	0%
	Manual	6	100%
<input checked="" type="checkbox"/>	DERS	6	100%
	Basic Mode	0	0%

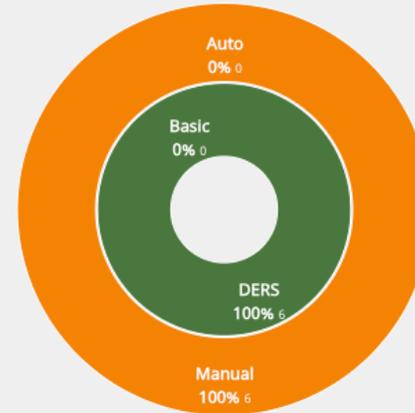

Care Area	Programming Mode	Dose Change
<input checked="" type="checkbox"/> All Selected Care Areas	All	6
<input type="checkbox"/> Adult ICU	All	6

Reports

1. Arrow should change direction when expanded
2. Line should stretch to width if equation

Quick Search

Search By

01-JAN-2023 - 02-OCT-2023

APPLY

Start Date: 01-JAN-2023

End Date: 02-OCT-2023

[Clear Filters](#)

Advanced Filters ▾

[Save / Manage Filters](#)

Compliance

Limits

Dose Rate Change

Device Usage

[← BACK TO DOSE RATE CHANGE REPORT](#)

All Selected Drugs

All Selected Care Areas

Date / Time	Drug Name (Modifier)	Concentration	Program Type	Mode Type	Change Requested	Result	Care Area
19-JUL-2023, 00:50:20	ceFAZolin	1 g / 50 mL	Manual	DERS	+172.7%	Within Limits	Adult ICU
19-JUL-2023, 00:50:03	ceFAZolin	1 g / 50 mL	Manual	DERS	-60%	Within Limits	Adult ICU
19-JUL-2023, 00:28:17	ceFAZolin	1 g / 50 mL	Manual	DERS	+66.7%	Within Limits	Adult ICU
19-JUL-2023, 00:25:01	ceFAZolin	1 g / 50 mL	Manual	DERS	-26.7%	Within Limits	Adult ICU
19-JUL-2023, 00:11:04	ceFAZolin	1 g / 50 mL	Manual	DERS	-55%	Within Limits	Adult ICU
19-JUL-2023, 00:09:10	ceFAZolin	1 g / 50 mL	Manual	DERS	+10%	Within Limits	Adult ICU

10 ▾

Items Per Page

< 1 >

6 Items

EXPORT

Reports

1. Green hard to read

< BACK TO DOSE RATE CHANGE REPORT

All Selected Drugs

All Selected Care Areas

Date / Time ^	Drug Name (Modifier)	Concentration	Program Type	Mode Type	Change Requested	Result	Care Area										
19-JUL-2023, 00:09:10	ceFAZolin	1 g / 50 mL	Manual	DERS	+10%	Within Limits	Adult ICU										
<div><p>HOSPITAL DEVICE ID 1 PUMP SERIAL NUMBER 02BF899A8BD6-1</p><p>LOCATION / ACCESS POINT DEVICE MODEL & TYPE 3 Novum IQ LVP</p><p>INFUSION ID 000000001689747345297204600000 ENTERPRISE GROUP NPP</p></div> <div><p>Dose Rate Change Parameters 2 View Infusion Story</p><table border="1"><thead><tr><th>DOSE RATE-PRIOR</th><th>DOSE RATE-AFTER</th><th>CHANGE</th><th>CHANGE LIMITS</th><th>DATE / TIME</th></tr></thead><tbody><tr><td>90.909 mL/hr</td><td>100 mL/hr</td><td>+10%</td><td>% / %</td><td>19-JUL-2023, 00:09:10</td></tr></tbody></table></div>								DOSE RATE-PRIOR	DOSE RATE-AFTER	CHANGE	CHANGE LIMITS	DATE / TIME	90.909 mL/hr	100 mL/hr	+10%	% / %	19-JUL-2023, 00:09:10
DOSE RATE-PRIOR	DOSE RATE-AFTER	CHANGE	CHANGE LIMITS	DATE / TIME													
90.909 mL/hr	100 mL/hr	+10%	% / %	19-JUL-2023, 00:09:10													
19-JUL-2023, 00:11:04	ceFAZolin	1 g / 50 mL	Manual	DERS	-55%	Within Limits	Adult ICU										
19-JUL-2023, 00:25:01	ceFAZolin	1 g / 50 mL	Manual	DERS	-26.7%	Within Limits	Adult ICU										
19-JUL-2023, 00:28:17	ceFAZolin	1 g / 50 mL	Manual	DERS	+66.7%	Within Limits	Adult ICU										
19-JUL-2023, 00:50:03	ceFAZolin	1 g / 50 mL	Manual	DERS	-60%	Within Limits	Adult ICU										
19-JUL-2023, 00:50:20	ceFAZolin	1 g / 50 mL	Manual	DERS	+172.7%	Within Limits	Adult ICU										

Reports

1. Background color should be same as row
2. Outline is difficult to see
3. Too much vertical space between rows

Library

NovumLVPSYRUS_Demo

[Change Drug Library](#)

Quick Search

Search By

[01-JAN-2023 - 02-OCT-2023](#)

APPLY

Start Date: 01-JAN-2023

End Date: 02-OCT-2023

[Clear Filters](#)

Advanced Filters

[Save / Manage Filters](#)

Compliance Limits Dose Rate Change **Device Usage**

View : CONNECTIVITY UTILIZATION

TIME OF DAY DAY OF WEEK SEASONAL

Serial Number	Connection %	Device Status	Avg. Utilization	Infusions	Errors	Last Updated	Certificate Exp.
No data available							
<input type="text" value="10"/>	Items per page	1	<input type="button" value="1"/>	0 items	<input type="button" value="EXPORT"/>		

Reports

1. Should be inline

Can't judge too much with no data

Library

NovumLVPSYRUS_Demo

[Change Drug Library](#)

Quick Search

Search By

[01-JAN-2023 - 02-OCT-2023](#)

APPLY

Start Date: 01-JAN-2023

End Date: 02-OCT-2023

[Clear Filters](#)

Advanced Filters

[Save / Manage Filters](#)

Compliance

Limits

Dose Rate Change

Device Usage

View :

CONNECTIVITY

UTILIZATION

TIME OF DAY

DAY OF WEEK

SEASONAL

Serial Number	Connection %	Device Status	Avg. Utilization	Infusions	Errors	Last Updated	Certificate Exp.
---------------	--------------	---------------	------------------	-----------	--------	--------------	------------------

No data available

10

Items per page

1



1



0 items

EXPORT

Reports

1. Should be inline

Can't judge too much with no data

Patient Index Manager

You are not authorized to view this page!

Patient Index Manager

If login does not allow it, page should not be in menu or should at least be grayed out